

Protein Calibration Curves on BioSep-SEC-S (2)

Column: BioSep™ 5 µm SEC-s3000 290 Å, LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-2146-K0

Elution Type: Isocratic

Eluent A: 100mM Sodium Phosphate, 300mM NaCl, pH 7.0

Gradient

Step No.	Time (min)	Pct A
1	0	100


Flow Rate: 1 mL/min

Col. Temp.: ambient



Detection: Laser Polarimeter @ 0.0000000000 (ambient)

Analyst Note: Application Focus: To show the calibration curves for the three different BioSep-SEC phases

For this application, calibration curves are generated for the three different BioSep-SEC phases using standard proteins of defined molecular weights and non-denaturing operating conditions (100 mM phosphate, 300 mM NaCl, pH 7.0, 1 mL/min). For each media (300 x 7.8 mm x 2 µm). Specific applications used to generate curves were: App ID# 18892 for BioSep2000, App ID# 18928 for BioSep3000 and App ID# 18932 for BioSep4000. While a polynomial was used to generate calibration curves in this example, one could use a linear plot if one removed the data points

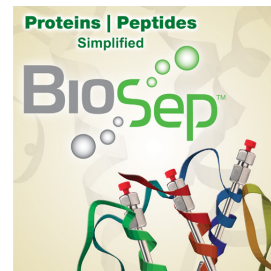

Products used in this application:

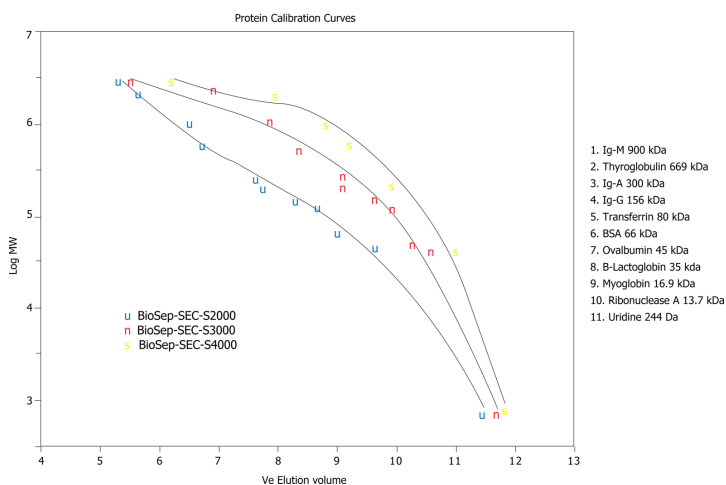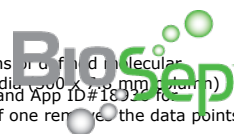

ANALYTES:

1 11 Proteins

