

## Linear MW Calibration on BioSep-SEC-S3000(2), Med Range MW

**Column:** BioSep™ 5 µm SEC-s3000 290 Å, LC Column 300 x 7.8 mm, Ea

**Dimensions:** 300 x 7.8 mm ID

**Order No:** 00H-2146-K0

**Elution Type:** Isocratic

**Eluent A:** 10mM Tris, 150mM NaCl, pH 7.4 □□

| Gradient | Step No. | Time (min) | Pct A |
|----------|----------|------------|-------|
| Profile: | 1        | 0          | 100   |

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** Laser Polarimeter @ 0.0000000000 (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, GFC 3000 4 x 3.0mm ID, 10/Pk Part No.: AJ0-4488

- Holder Part No.: KJ0-4282



Products used in this application:

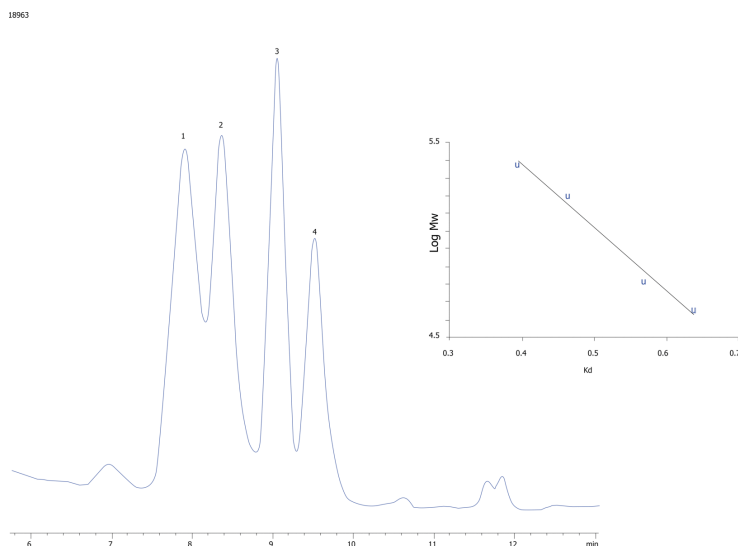

### ANALYTES:

- 1** IgA  
Retention Time: 7.9 min
- 2** beta-Amylase  
Retention Time: 8.36 min
- 3** BSA  
Retention Time: 9.05 min
- 4** Ovalbumin  
Retention Time: 9.51 min

