

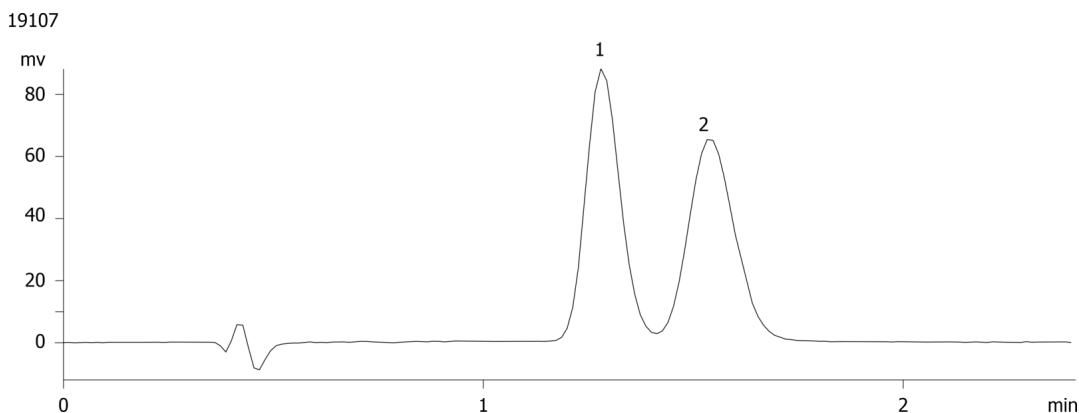
Impact of column length on Chiral Separations: 50mm length**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Ethanol with 0.1% DEA**Eluent B:** Carbon Dioxide

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	25	75

Flow Rate: 2.5 mL/min**Col. Temp.:** 35 °C**Detection:** Laser Polarimeter @ 0.0000000000 (35 °C)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Silica 4 x 3.0mm, 10/Pk Part No.: AJ0-4348

- Holder Part No.: KJ0-4282

**ANALYTES:****1** Propafenone