


BSA tryptic map in formic acid using Kinetex 2.6u C18 150x4.6mm**Column:** Kinetex® 2.6 µm C18 100 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4462-E0**Elution Type:** Gradient**Eluent A:** Water with 0.1% formic acid**Eluent B:** Acetonitrile with 0.1% formic acid

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	97	3
	2	30	50	50

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)**Products used in this application:****ANALYTES:**

- 1 Bovine Serum Albumin (BSA)

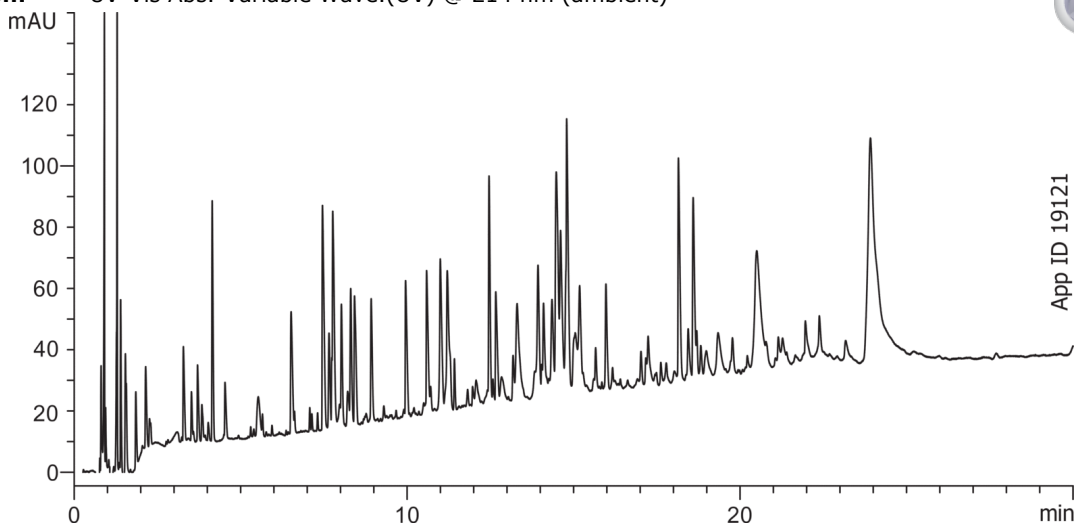