

CD Detection of TSO on LUX 5µm Cellulose-1 250 x 4.6 mm ID in NP

Column: Lux® 5 µm Cellulose-1, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4459-E0

Elution Type: Isocratic

Eluent A: 0.1% DEA in Hexane

Eluent B: 0.1% DEA in IPA

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	10	90	10

Flow Rate: 1 mL/min

Col. Temp.: 25 °C


Detection: Circular Dichroism @ 0.0000000000 (25 °C)

Analyst Note: Please see the corresponding CD trace (attached document).

19149



Products used in this application:

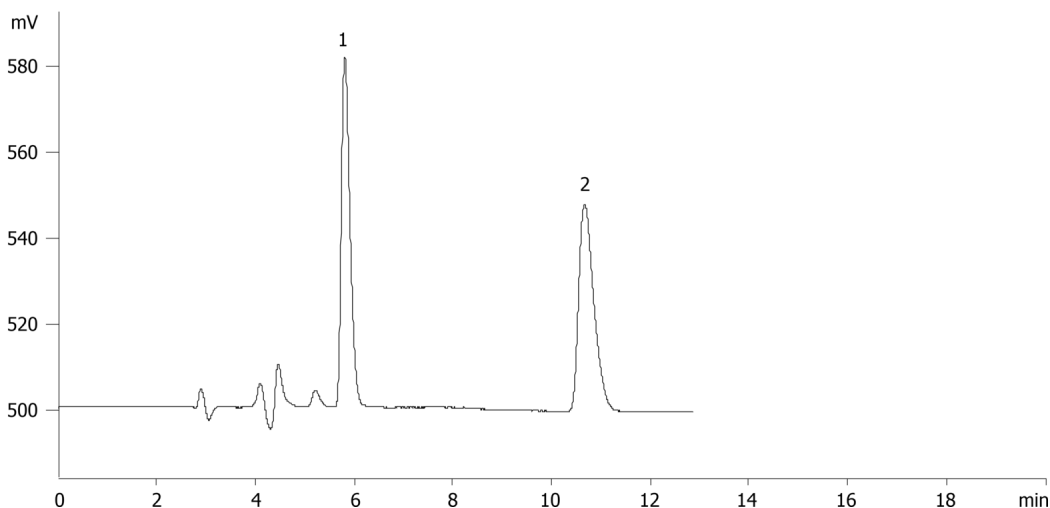

ANALYTES:

- 1 trans-Stilbene oxide
- 2 trans-Stilbene oxide

