

Chiral HPLC Purification of Atenolol on Lux Cellulose-1

Column: Lux® 5 µm Cellulose-1, LC Column 250 x 21.2 mm, AXIA™ Packed, Ea

Dimensions: 250 x 21.2 mm ID

Order No: 00G-4459-P0-AX

Elution Type: Isocratic

Eluent A: 0.1% DEA in Ethanol

Eluent B: 0.1% DEA in Hexane

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	1	10	20	80

Flow Rate: 50 mL/min

Col. Temp.: 20 °C

Detection: Circular Dichroism @ 0.0000000000 (20 °C)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard PREP Cartridge, Lux Cellulose-1, 15 x 21.2mm ID, Ea Part No.: AJ0-8405

- Holder Part No.: AJ0-8223

19269



Products used in this application:

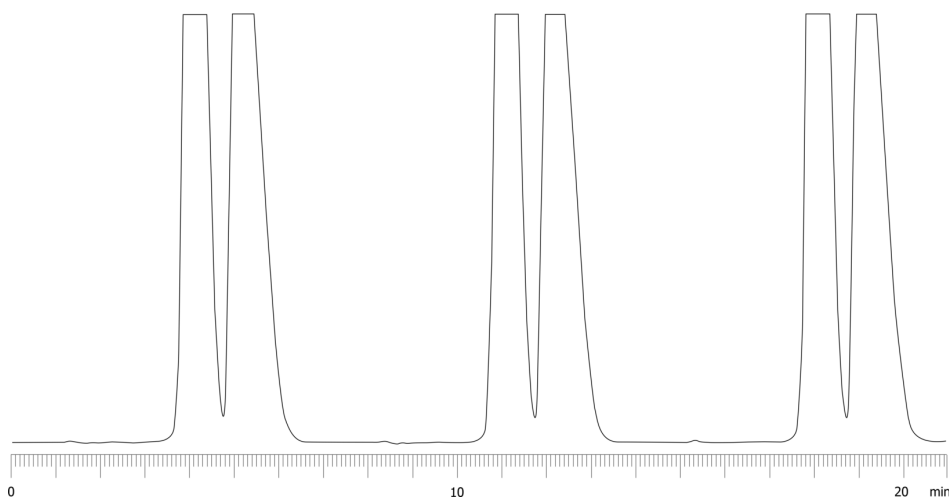

ANALYTES:

1 Atenolol

Retention Time: 4 min

2 Atenolol

Retention Time: 5.5 min

