


**Organophosphorous Pesticides on ZB-MR-1 by MSD**

**Column:** Zebron™ ZB-MultiResidue™-1, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea  
**Phase:** Proprietary Pesticides Phase  
**Dimensions:** 30 meters x 0.25 mm x 0.25 µm  
**Order No:** 7HG-G016-11  
**Oven Profile:**

**Carrier Gas:** Constant Flow Helium, 0.9 mL/min  
**Injection:** Splitless 0.000000000:1 0 µL @ 260°C  
**Detection:** Refractive Index (320°C)  
**Analyst Note:** Application Courtesy of Dr. Bernard Wronski of Landeslabor Brandenburg.



Products used in this application:

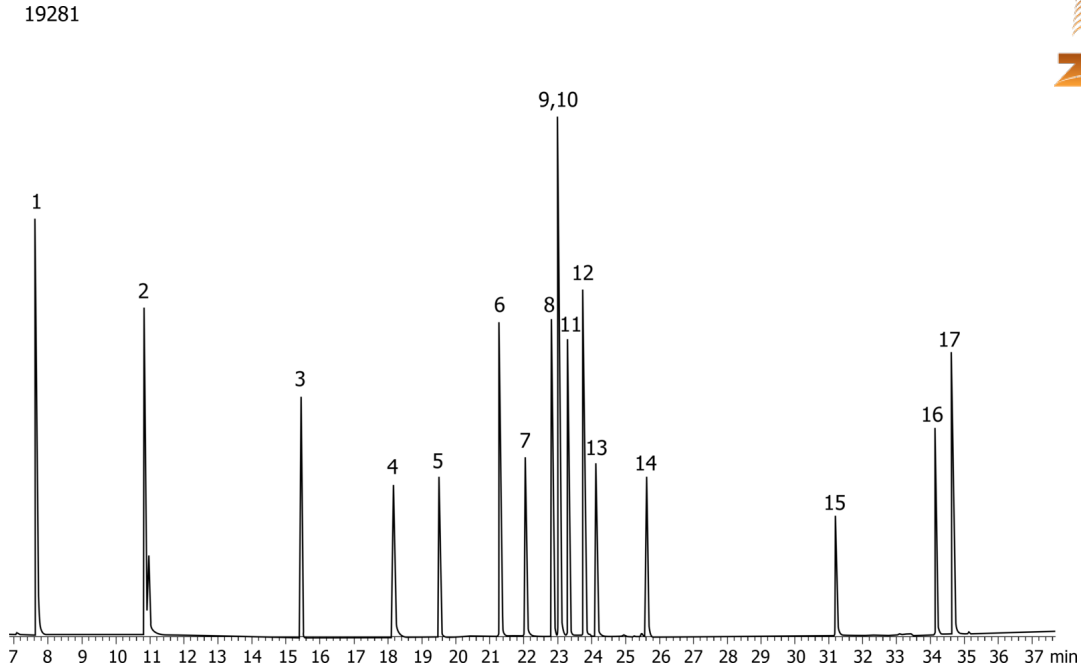

## Organophosphorous Pesticides on ZB-MR-1 by MSD

### ANALYTES:

- 1 Dichlorvos
- 2 Mevinphos
- 3 Trifluralin
- 4 Dimethoate
- 5 Etrimfos
- 6 Chlorpyrifos-methyl
- 7 Parathion-methyl
- 8 Malathion-d10
- 9 Malathion
- 10 Fenitrothion
- 11 Chlorpyrifos-ethyl
- 12 Fenthion
- 13 Parathion
- 14 Chlorfenvinphos
- 15 Triazophos
- 16 Azinphos-methyl
- 17 Azinphos-ethyl

