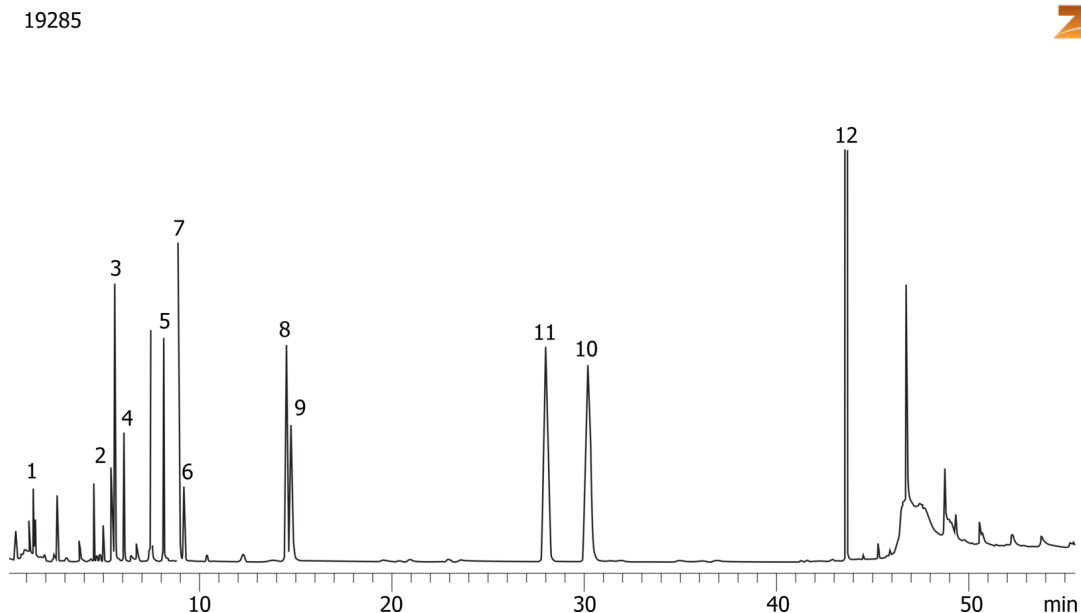


EPA Method 552.2 - Haloacetic acids on a ZB-XLB-HT (Dual column 0.32 ID)**Column:** Zebron™ ZB-XLB-HT, GC Cap. Column 30 m x 0.32 mm x 0.25 µm, Ea**Phase:** Proprietary XLB Phase**Dimensions:** 30 meters x 0.32 mm x 0.25 µm**Order No:** 7HM-G024-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 2.5 mL/min**Injection:** Split 25.0000000000:1 0 µL @ 250°C**Detection:** Refractive Index (350°C)**Analyst Note:** Analytes are methylated and 100 ppb in isooctane.

Dual column with ZB-35-HT

**Products used in this application:**

EPA Method 552.2 - Haloacetic acids on a ZB-XLB-HT(Dual column 0.32 ID)

ANALYTES:

- 1 Chloroacetic acid
- 2 Bromoacetic acid
- 3 Dichloroacetic acid
- 4 Dalapon
- 5 Trichloroacetic acid
- 6 1,2,3-Trichloropropane
- 7 Bromochloroacetic acid
- 8 Bromodichloroacetic acid
- 9 Dibromoacetic acid
- 10 Chlorodibromoacetic acid
- 11 2,3-Dibromopropionate
- 12 Tribromoacetic acid

