


USP Glycerin Purity Assay on ZB-BAC-2**Column:** Zebtron™ ZB-BAC-2, GC Cap. Column 30 m x 0.53 mm x 2.00 µm, Ea**Phase:** Proprietary Blood Alcohol Phase**Dimensions:** 30 meters x 0.53 mm x 2 µm**Order No:** 7HK-G022-32**Oven Profile:****Carrier Gas:** Constant Flow Helium, 4.5 mL/min**Injection:** Split 10.0000000000:1 0 µL @ 220°C**Detection:** Refractive Index (250°C)**Analyst Note:** Concentrations were 2.0 mg/mL Glycerin, 0.050 mg/mL ethylene glycol, 0.050 mg/mL diethylene glycol, and 0.10 mg/mL 2,2,2-trichloroethanol (IS)
Run as a dual column analysis with ZB-BAC-1
Resolution = 11.09 between glycerin and diethylene glycol**Products used in this application:**

19291

**ANALYTES:**

- 1 Ethylene glycol
- 2 2,2,2-Trichloroethanol
- 3 Diethylene glycol
- 4 Glycerin

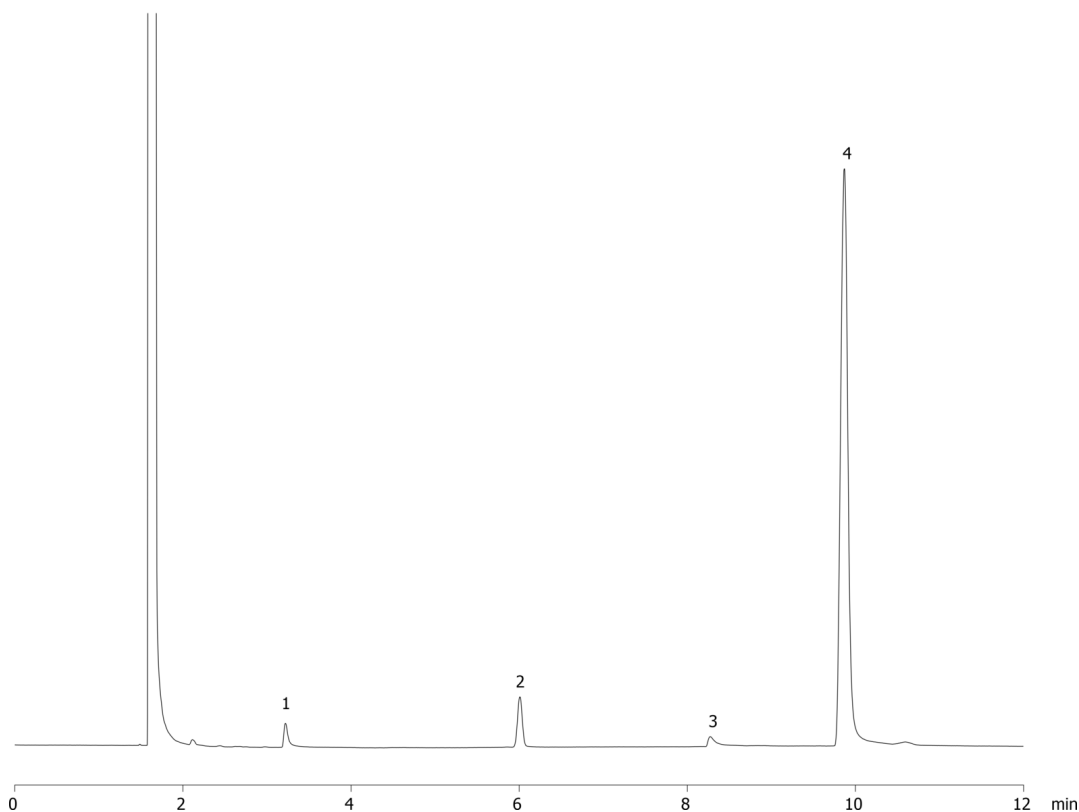