


Valproic Acid Assay (USP) on ZB-FFAP, 30 m x 0.25 mm x 0.25 µm**Column:** Zebron™ ZB-FFAP, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Nitroterephthalic Acid Modified Polyethylene Glycol**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G009-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 3 mL/min**Injection:** Split 30.0000000000:1 0 µL @ 200°C**Detection:** Refractive Index (250°C)**Analyst Note:** Valproic acid was 1.0 mg/mL and Nonanoic acid was 1.2 mg/mL both in hexane.**Products used in this application:****ANALYTES:**

- 1 Valproic acid
- 2 Nonanoic acid

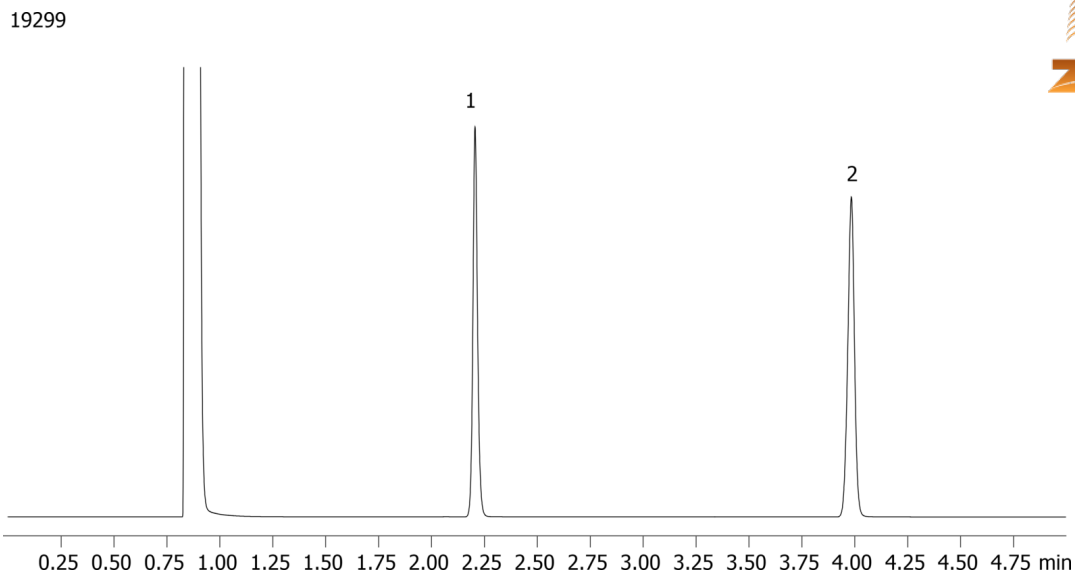