


Valproic Acid Impurities (USP) on ZB-FFAP, 30 m x 0.32 mm x 0.25 µm**Column:** Zebron[™] ZB-FFAP, GC Cap. Column 30 m x 0.32 mm x 0.25 µm, Ea**Phase:** Nitroterephthalic Acid Modified Polyethylene Glycol**Dimensions:** 30 meters x 0.32 mm x 0.25 µm**Order No:** 7HM-G009-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.1 mL/min**Injection:** Split 100.000000000:1 0 µL @ 240°C**Detection:** Refractive Index (260°C)**Analyst Note:** Concentrations were 1.0 mg/mL butyric acid, 1.0 mg/mL valeric acid, 0.1 mg/mL related compound A in Valproic**Products used in this application:****ANALYTES:**

- 1 Butyric acid
- 2 Valeric acid
- 3 Valproic acid
- 4 Valproic acid related compound A

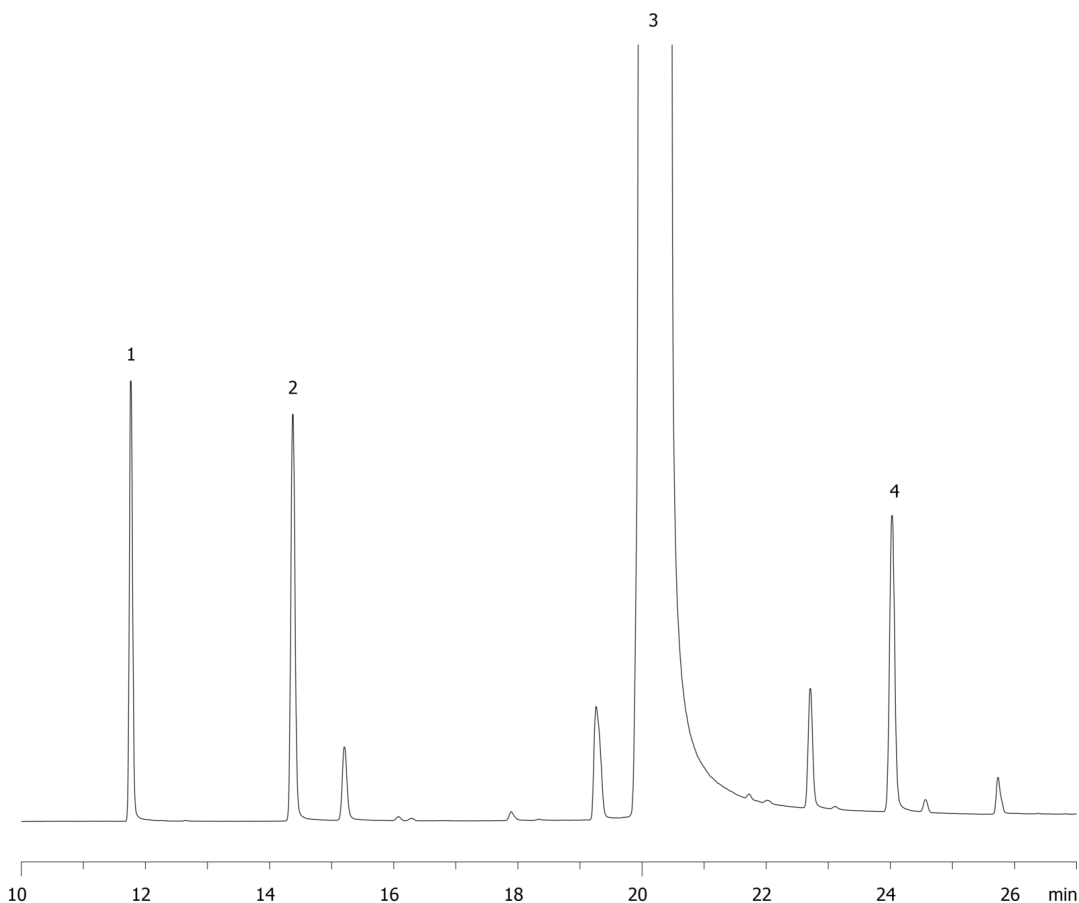