

Steroids Mix using Kinetex 2.6µm C18 50x2.1mm

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4462-AN

Elution Type: Gradient

Eluent A: water

Eluent B: ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	6	40	60
	3	6.01	80	20
	4	8	80	20


Flow Rate: 0.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)



Products used in this application:

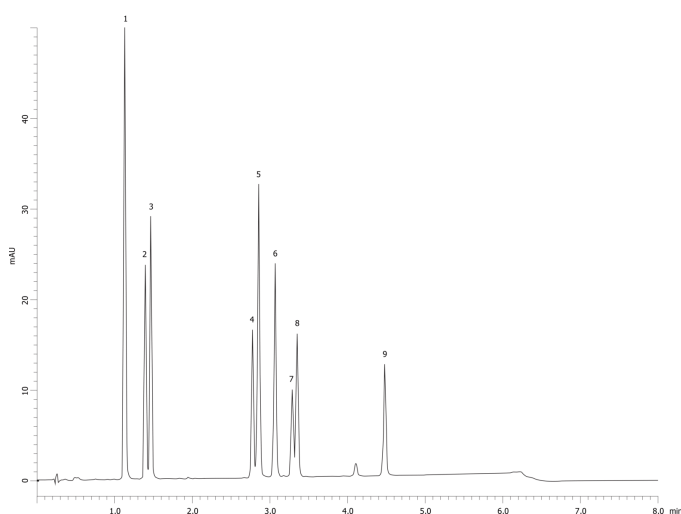

ANALYTES:

- 1 Estriol
Retention Time: 1.13 min
- 2 Hydrocortisone
Retention Time: 1.394 min
- 3 Cortisone
Retention Time: 1.464 min
- 4 Cortisone-21-acetate
Retention Time: 2.778 min
- 5 Estradiol
Retention Time: 2.855 min
- 6 21-Hydroxyprogesterone
Retention Time: 3.068 min
- 7 Estrone
Retention Time: 3.29 min
- 8 17-Hydroxyprogesterone
Retention Time: 3.352 min
- 9 Deoxycorticosterone acetate
Retention Time: 4.48 min

