

Amphetamines by LC/MS/MS using a Kinetex 2.6µm C18 50x2.1mm

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4462-AN

Elution Type: Gradient

Eluent A: 5mM Ammonium formate with 0.1% FA

Eluent B: Methanol with 0.1% FA

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	1	30	70
	3	3	30	70
	4	3	30	70

Flow Rate: 400

Col. Temp.: ambient

Detection: Circular Dichroism @ 0.0000000000 (ambient)

Analyst Note: SecurityGuard™ ULTRA Guard Cartridge System extends column lifetime.


- SecurityGuard ULTRA Cartridges, UHPLC C18 for 2.1mm ID Columns, 3/Pk Part No.: AJ0-8782

- Holder Part No.: AJ0-9000

19490



Products used in this application:

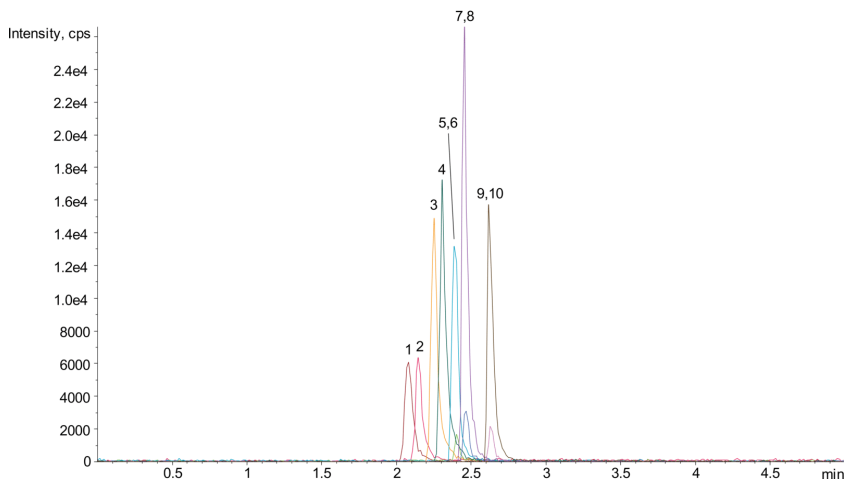

Amphetamines by LC/MS/MS using a Kinetex 2.6µm C18 50x2.1mm

ANALYTES:

- 1** Amphetamine-d11
- 2** Amphetamine
- 3** Methamphetamine-d14
- 4** Methamphetamine
- 5** MDA-d5
- 6** MDA
- 7** MDMA-d5
- 8** MDMA
- 9** MDEA-d5
- 10** MDEA

