


Active Drug Mix on ZB-5MS 30 x 0.25 x 0.25**Column:** Zebtron[™] ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G010-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 40 cm/sec**Injection:** On-Column 0.000000000:1 0 µL @ 163°C**Detection:** Refractive Index (320°C)**Analyst Note:** Active drug mix, 100 ppm in methanol

19582

**Products used in this application:****ANALYTES:**

- 1 Ibuprofen
- 2 Glutethimide
- 3 Dimenhydrinate
- 4 Desipramine
- 5 Oxymorphone
- 6 Quinidine
- 7 Alprazolam

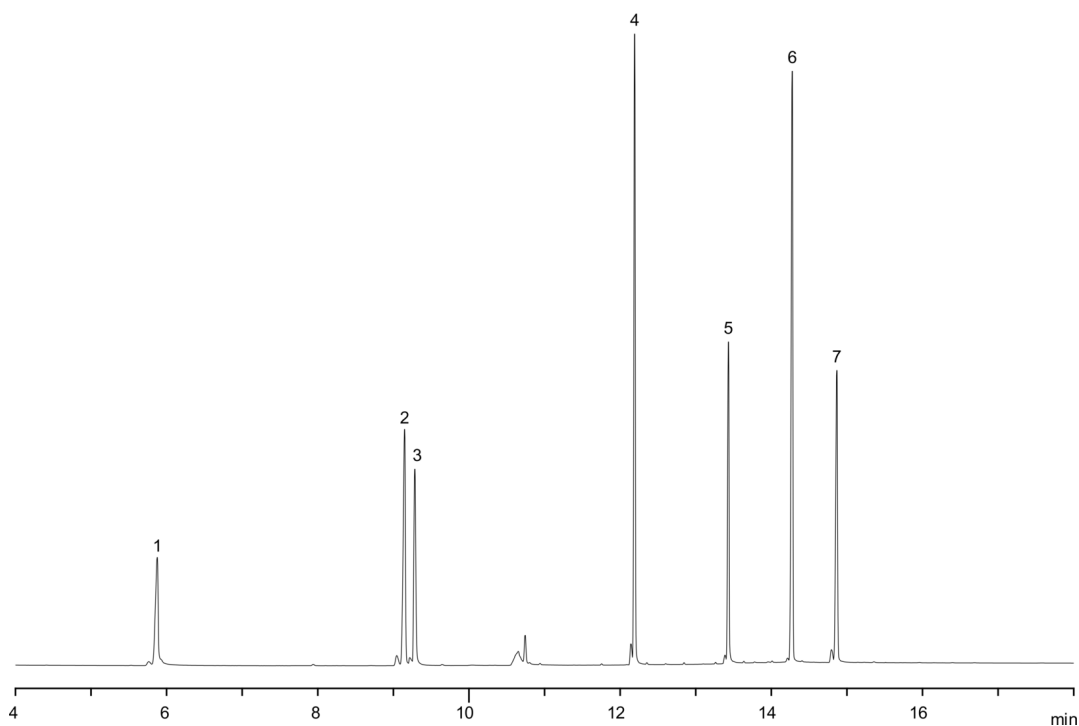