

## Ultra Inert Mix on ZB-5MS, 30 x 0.25 x 0.25

**Column:** Zebtron™ ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea

**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane

**Dimensions:** 30 meters x 0.25 mm x 0.25 µm

**Order No:** 7HG-G010-11

**Oven Profile:**

**Carrier Gas:** Constant Flow Helium, 40 cm/sec


**Injection:** On-Column 0.000000000:1 0 µL @ 68°C

**Detection:** Refractive Index (320°C)

**Analyst Note:** Ultra Inert mix, 100 ppm in acetone



Products used in this application:

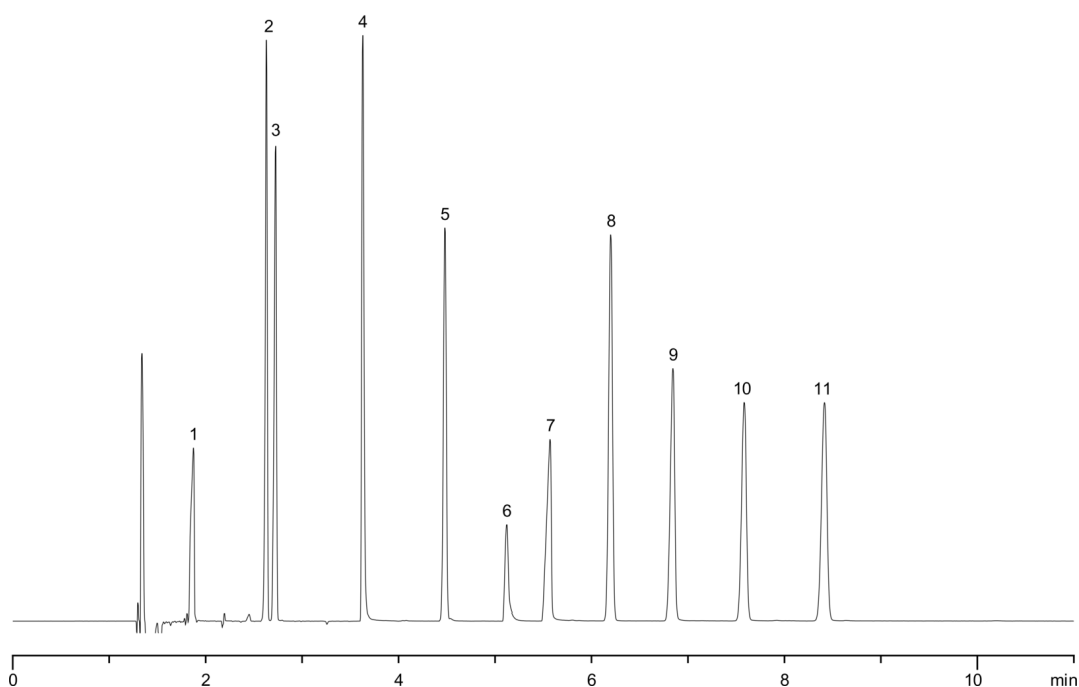

## Ultra Inert Mix on ZB-5MS, 30 x 0.25 x 0.25

### ANALYTES:

- 1** 1-Propionic acid
- 2** 1-Octene
- 3** n-Octane
- 4** 4-Methylpyrimidine
- 5** n-Nonane
- 6** Trimethyl phosphate
- 7** 1,2-Pentanediol
- 8** N-propylbenzene
- 9** 1-Heptanol
- 10** 3-Octanone
- 11** n-Decane

