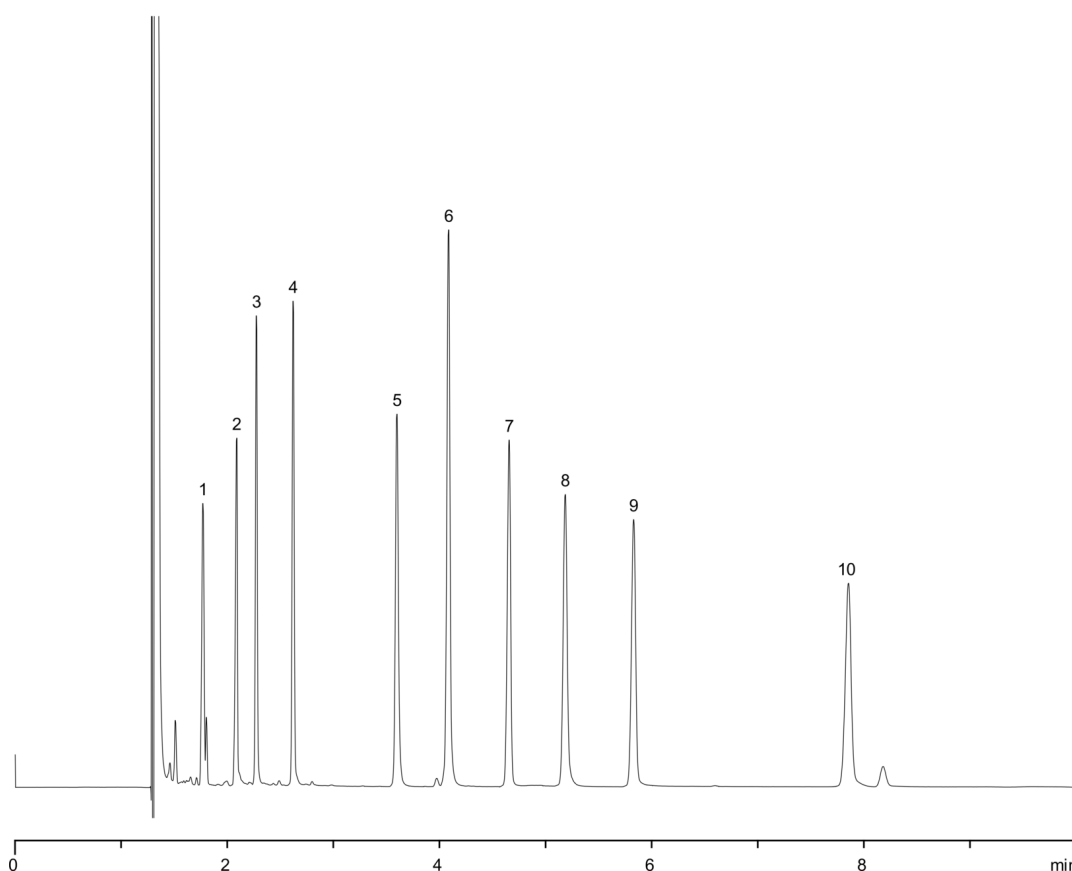


ZB-5MS Test Mix on ZB-5MS 30 x 0.25 x 0.25**Column:** Zebtron[™] ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G010-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 40 cm/sec**Injection:** On-Column 0.000000000:1 0 µL @ 143°C**Detection:** Refractive Index (320°C)**Analyst Note:** Analytes dissolved at 100 ppm in acetone

19587

**Products used in this application:**

ZB-5MS Test Mix on ZB-5MS 30 x 0.25 x 0.25

ANALYTES:

- 1 Decane
- 2 2-Ethylhexanoic acid
- 3 1,6-Hexanediol
- 4 4-Chlorophenol
- 5 Tridecane
- 6 1-Methylnaphthalene
- 7 1-Undecanol
- 8 Tetradecane
- 9 Dicyclohexylamine
- 10 Pentadecane

