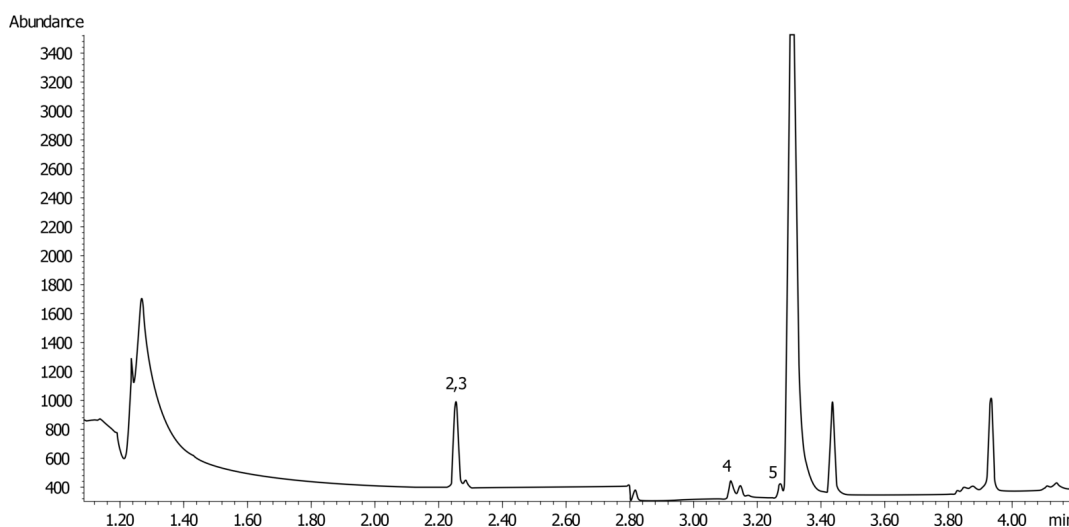


Dispersants in seawater on a ZB-Waxplus 30x0.25x0.25 by GCMS**Column:** Zebtron™ ZB-WAXplus™, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G013-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.5 mL/min**Injection:** Split 10.0000000000:1 0 µL @ 250°C**Detection:** Refractive Index (250°C)**Analyst Note:** Analytes are 1.7ppm in seawater

Mass spectrometer is operated in SIM mode: Window 1 (87 (2-butoxyethanol), 89 (IS)), Window 2 (propyleneglycol and IS (31))

**Products used in this application:**

19604

**ANALYTES:**

- 1 2-Butoxyethanol-d4
- 2 2-Butoxyethanol
- 3 Propylene glycol
- 4 Ethylene glycol (EG)

