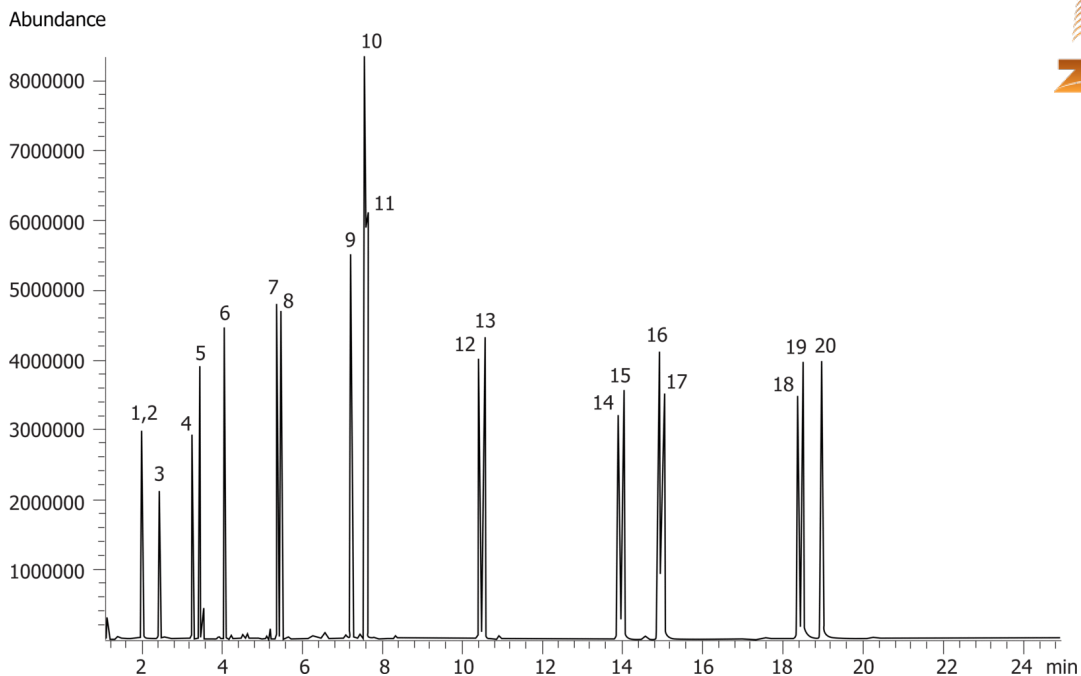


**Polycyclic Aromatic Hydrocarbons (PAHs) in Water by Zebron ZB-5MS****Column:** Zebron<sup>™</sup> ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25  $\mu$ m, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25  $\mu$ m**Order No:** 7HG-G010-11**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.4 mL/min**Injection:** Split 15.0000000000:1 0  $\mu$ L @ 310°C**Detection:** Refractive Index (270°C)

19651

**Products used in this application:**

**Polycyclic Aromatic Hydrocarbons (PAHs) in Water by Zebron ZB-5MS****ANALYTES:**

- 1 Naphthalene-d8
- 2 Naphthalene
- 3 2-Methylnaphthalene
- 4 Acenaphthalene
- 5 Acenaphthene
- 6 Fluorene
- 7 Anthracene
- 8 Phenanthrene
- 9 Fluoranthene
- 10 Pyrene-d10
- 11 Pyrene
- 12 Benzo[a]anthracene
- 13 Chrysene
- 14 Benzo[b]fluoranthene
- 15 Benzo[k]fluoranthene
- 16 Benzo[a]pyrene-d12
- 17 Benzo[a]pyrene
- 18 Indeno[1,2,3-cd]pyrene
- 19 Dibenzo[a,h]anthracene
- 20 Benzo[g,h,i]perylene

