

Protein Mixture Aeris 3.6µm XB-C18 Wide Pore 150x2.1mm

Column: Aeris™ 3.6 µm WIDEPORE XB-C18 200 Å, LC Column 150 x 2.1 mm, Ea

Dimensions: 150 x 2.1 mm ID

Order No: 00F-4482-AN

Elution Type: Gradient

Eluent A: 0.1% TFA in Water

Eluent B: 0.1% TFA in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	25	35	65

Flow Rate: 0.2 mL/min


Col. Temp.: 40 °C

Detection: Circular Dichroism @ 0.0000000000 (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC WIDEPORE C18 for 2.1mm ID, 3/Pk Part No.: AJ0-8783


- Holder Part No.: AJ0-9000

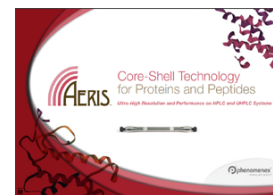

Products used in this application:



19771

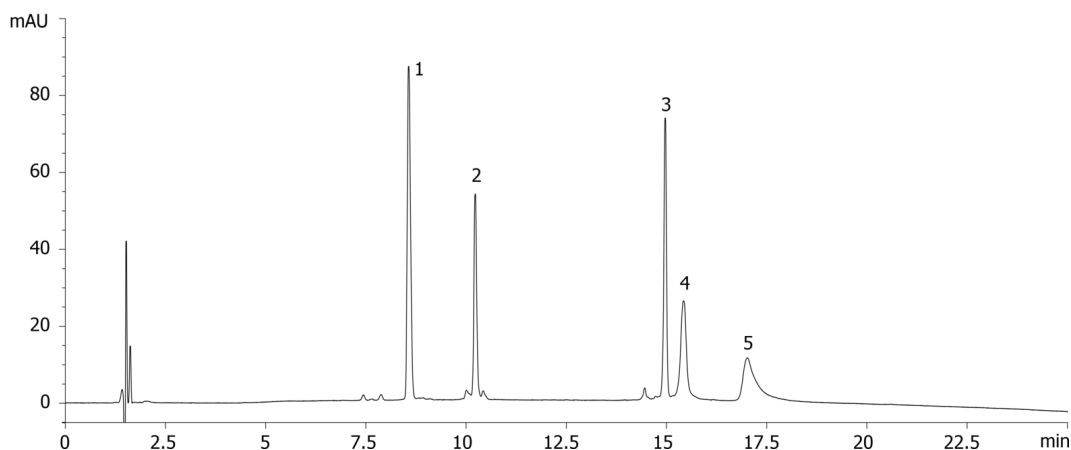

ANALYTES:

- 1** Angiotensin I
Retention Time: 8.565 min
- 2** Ribonuclease A
Retention Time: 10.226 min
- 3** Lysozyme
Retention Time: 14.967 min
- 4** Apotransferrin
Retention Time: 15.425 min
- 5** Bovine Serum Albumin (BSA)
Retention Time: 17.016 min

