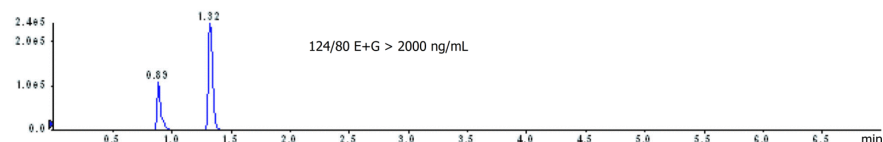
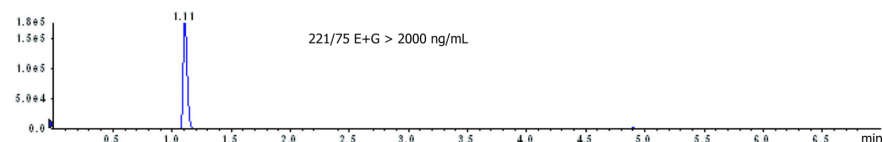
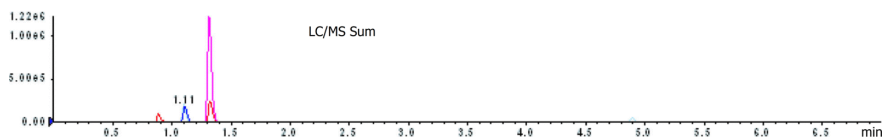


EtS EtG Sample #6 high levels**Column:** SynergiTM 2.5 μ m Hydro-RP 100 Å, LC Column 100 x 3 mm, Ea**Dimensions:** 100 x 3 mm ID**Order No:** 00D-4387-Y0**Elution Type:** Gradient**Eluent A:** 10mM Ammonium formate**Eluent B:** 50:50 Methanol: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	100	0
	2	4	20	80
	3	4.1	100	0
	4	7	100	0

Flow Rate: 0.5 mL/min**Col. Temp.:** ambient**Detection:** Mass Spectrometer (MS) @ amu (ambient)**Detector Info:** SCIEX<

19774

**ANALYTES:**

- 1 EtG (Ethyl glucuronide)
- 2 Ethyl glucuronide-d5 (EtG-d5)
- 3 EtS (Ethyl sulfate)
- 4 EtS (Ethyl sulfate)

