


Bioethanol using ZB-624 30x.25x1.40um**Column:** Zebron[™] ZB-624, GC Cap. Column 30 m x 0.25 mm x 1.40 µm, Ea**Phase:** 6% Cyanopropylphenyl 94% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 1.4 µm**Order No:** 7HG-G005-27**Oven Profile:****Carrier Gas:** Constant Flow Helium, 1.6 mL/min**Injection:** Split 50.0000000000:1 0 µL @ 140°C**Detection:** Refractive Index (200°C)**Analyst Note:** Solvent: Ethanol**Products used in this application:**

Bioethanol using ZB-624 30x.25x1.40um

ANALYTES:

- 1 Acetaldehyde
- 2 Methanol
- 3 Ethanol
- 4 Impurity from Ethanol
- 5 Propanol
- 6 Ethyl acetate
- 7 2-Butanol
- 8 2-Methylpropanol
- 9 Butanol
- 10 Acetal
- 11 3-Methylbutanol
- 12 2-Methylbutanol
- 13 4-Methylpentanol

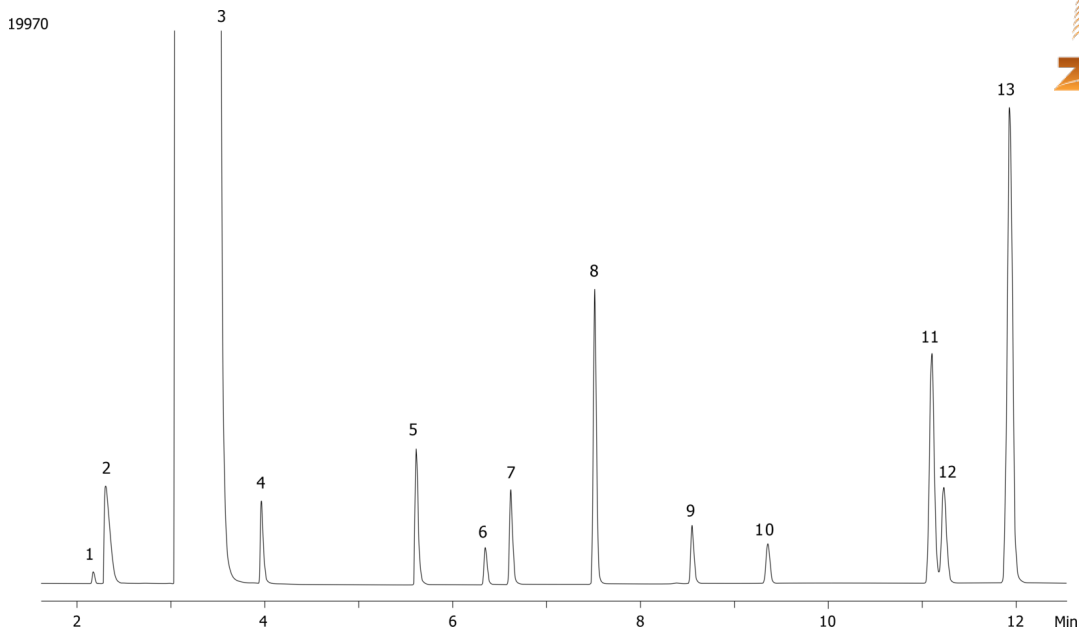