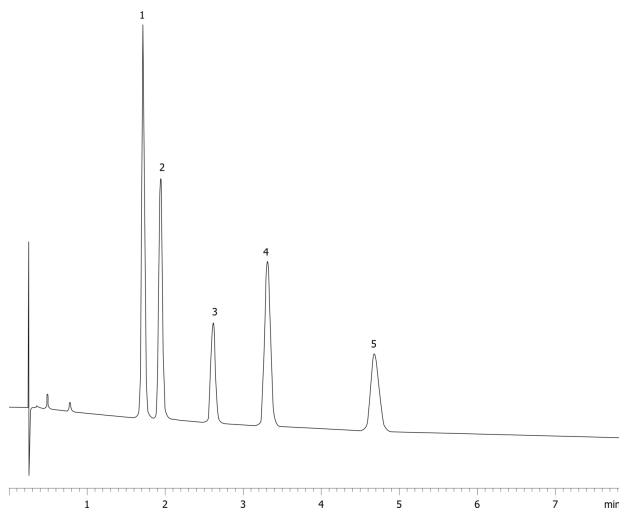


Substituted Methoxybenzenes on Kinetex 2.6 μ Phenyl-Hexyl**Column:** Kinetex® 2.6 μ m Phenyl-Hexyl 100 Å, LC Column 50 x 4.6 mm, Ea**Dimensions:** 50 x 4.6 mm ID**Order No:** 00B-4495-E0**Elution Type:** Isocratic**Eluent A:** Water**Eluent B:** MeOH

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	60	40

Flow Rate: 1.85 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Diode Array (PDA) @ 254 nm (ambient)

19985

**Products used in this application:****ANALYTES:**

- 1** 1,2-Dimethoxybenzene
Retention Time: 1.737 min
- 2** 1,2,4-Trimethoxybenzene
Retention Time: 1.963 min
- 3** 1,4-Dimethoxybenzene
Retention Time: 2.638 min
- 4** 1,3-Dimethoxybenzene
Retention Time: 3.34 min
- 5** 1,3,5-Trimethoxybenzene
Retention Time: 4.7 min

