

Separation of Benzodiazepines

Column: Synergi[™] 4 µm Hydro-RP 80 Å, LC Column 50 x 2 mm, Ea

Dimensions: 50 x 2 mm ID

Order No: 00B-4375-B0

Elution Type: Gradient

Eluent A: 0.1% Formic acid in DI H₂O

Eluent B: 0.1% formic Acid in ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0.5	90	10
	2	0.7	80	20
	3	5.9	25	75
	4	6	5	95
	5	6.2	5	95
	6	6.3	90	10
	7	8.5	90	10

Flow Rate: 500


Col. Temp.: 40 °C

Detection: Circular Dichroism @ 0.0000000000 (ambient)

Analyst Note: SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, AQ C18 4 x 2.0mm ID, 10/Pk Part No.: AJ0-7510

- Holder Part No.: KJ0-4282

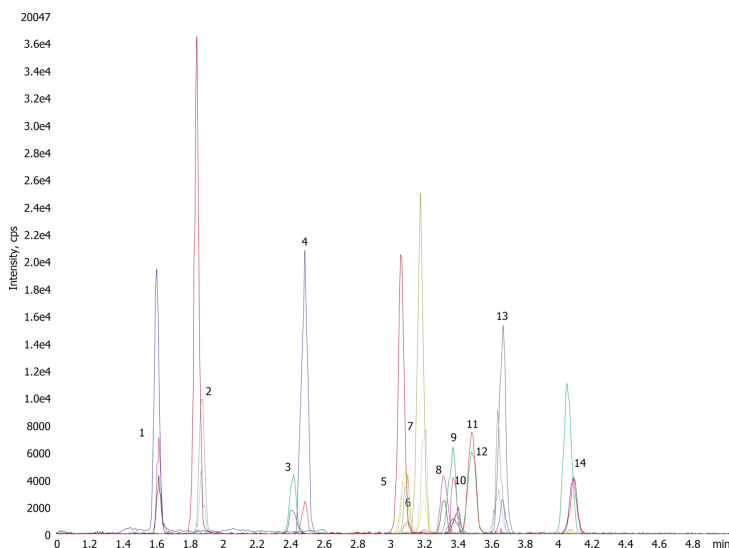

Separation of Benzodiazepines

ANALYTES:

- 1** 7-Aminoclonazepam
- 2** 7-Aminoflurazepam
- 3** Midazolam
- 4** Flurazepam
- 5** a-OH-Alprazolam
- 6** a-OH-Triazolam
- 7** Oxazepam
- 8** Lorazepam
- 9** Clonazepam
- 10** Nordiazepam
- 11** Alprazolam
- 12** Flunitrazepam
- 13** Temazepam
- 14** Diazepam

