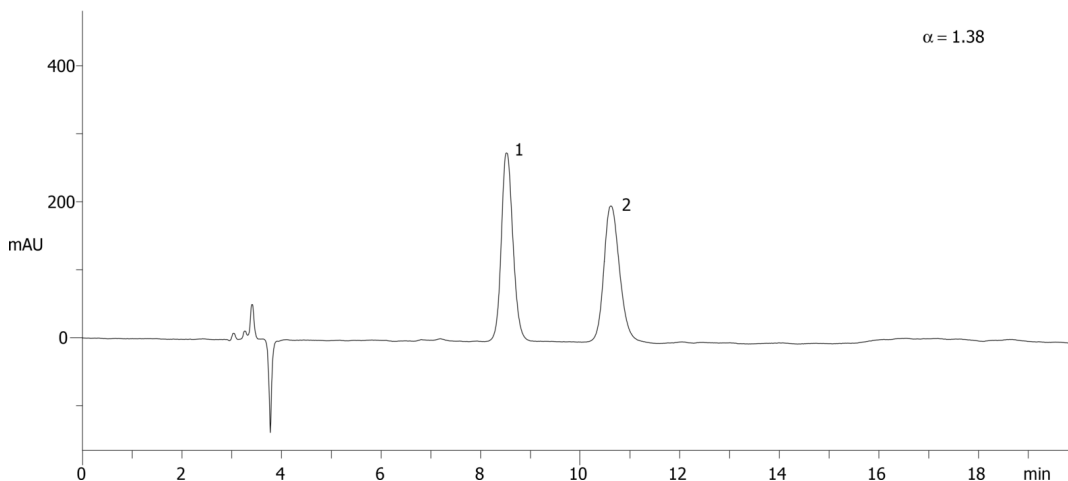


Ambucetamid on Lux 5 μ m Amylose-2 in NP**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Hexane / IPA (80:20) DEA (0.1%)**Gradient** **Step No.** **Time (min)** **Pct A****Profile:** **1** 0 100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** Circular Dichroism @ 0.0000000000 (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Amylose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8470

- Holder Part No.: KJ0-4282

20097

**ANALYTES:****1** Ambucetamid

Retention Time: 8.5 min

2 Ambucetamid

Retention Time: 10.6 min

