

Flumequine on Lux 5µm Cellulose-2 in RP

Column: Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4457-E0
Elution Type: Isocratic
Eluent A: ACN / 20mM NH₄HCO₃ (60:40) with DEA (0.1%) in buffer
Gradient


Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Circular Dichroism @ 0.0000000000 (ambient)
Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8366
- Holder Part No.: KJ0-4282



Products used in this application:

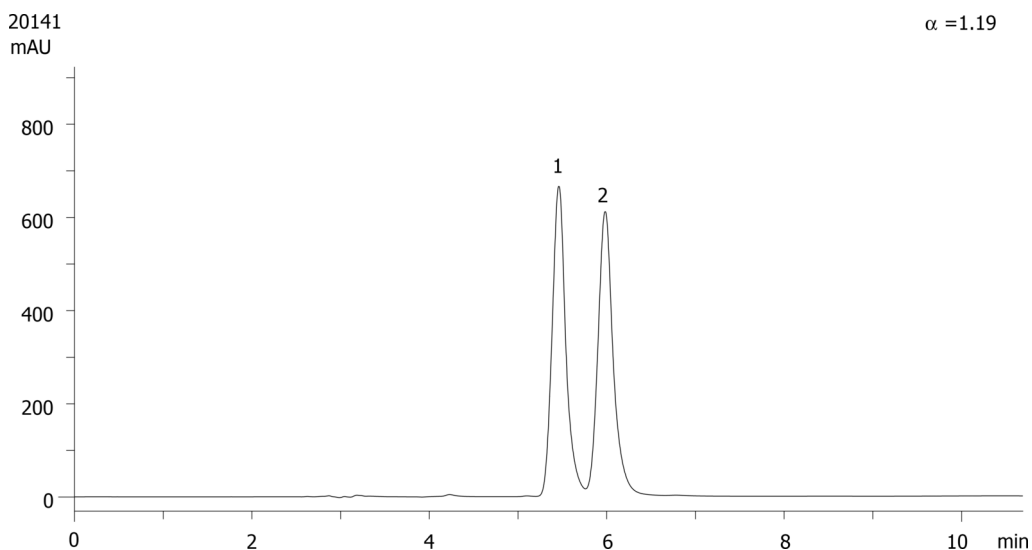

$\alpha = 1.19$

ANALYTES:

- 1** Flumequine
Retention Time: 5.46 min
- 2** Flumequine
Retention Time: 5.98 min

