

FMOC-Ile-OH on Lux 5µm Cellulose-4 in RP

Column: Lux® 5 µm Cellulose-4, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4491-E0
Elution Type: Isocratic
Eluent A: MeOH / 20mM NH₄HCO₃ (80:20) with Formic acid (0.1%) in buffer
Gradient


Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Circular Dichroism @ 0.0000000000 (ambient)
Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-4 4 x 3.0 mm, 10/Pk Part No.: AJ0-8627
 - Holder Part No.: KJ0-4282



Products used in this application:

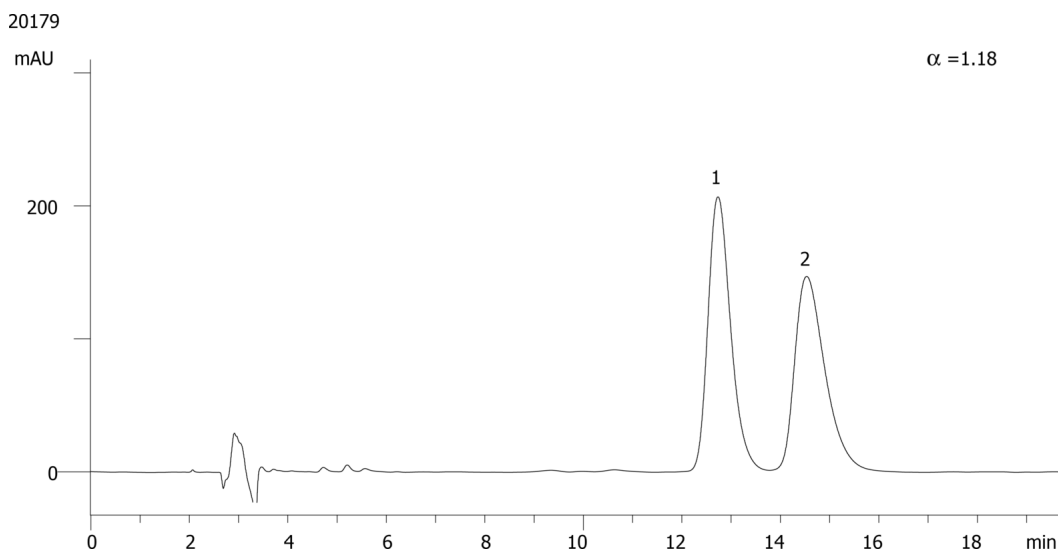

ANALYTES:

- 1 FMOC-Ile-OH
Retention Time: 12.8 min
- 2 FMOC-Ile-OH
Retention Time: 14.6 min

