

Physostigmine (Eserine) on Lux 5µm Cellulose-2 in RP

Column: Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4457-E0
Elution Type: Isocratic
Eluent A: MeOH / 20mM NH₄HCO₃+0.1%DEA (80:20) with None in buffer
Gradient


Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Circular Dichroism @ 0.0000000000 (ambient)
Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8366
 - Holder Part No.: KJ0-4282



Products used in this application:

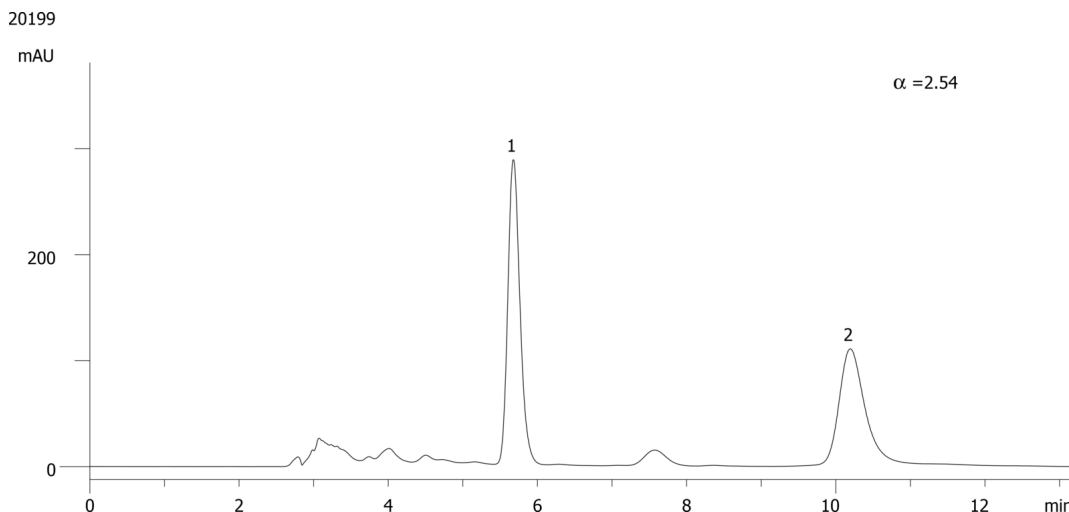

ANALYTES:

- 1 Physostigmine
Retention Time: 5.68 min
- 2 Physostigmine
Retention Time: 10.19 min

