

## Tropicamide on Lux 5µm Cellulose-2 in RP

**Column:** Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea  
**Dimensions:** 250 x 4.6 mm ID  
**Order No:** 00G-4457-E0  
**Elution Type:** Isocratic  
**Eluent A:** MeOH / 20mM NH<sub>4</sub>HCO<sub>3</sub>+0.1%DEA (60:40) with None in buffer  
**Gradient**

| Step No. | Time (min) | Pct A |
|----------|------------|-------|
| 1        | 0          | 100   |


**Profile:**  
**Flow Rate:** 1 mL/min  
**Col. Temp.:** ambient  
**Detection:** Circular Dichroism @ 0.0000000000 (ambient)  
**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8366
- Holder Part No.: KJ0-4282

20241



Products used in this application:

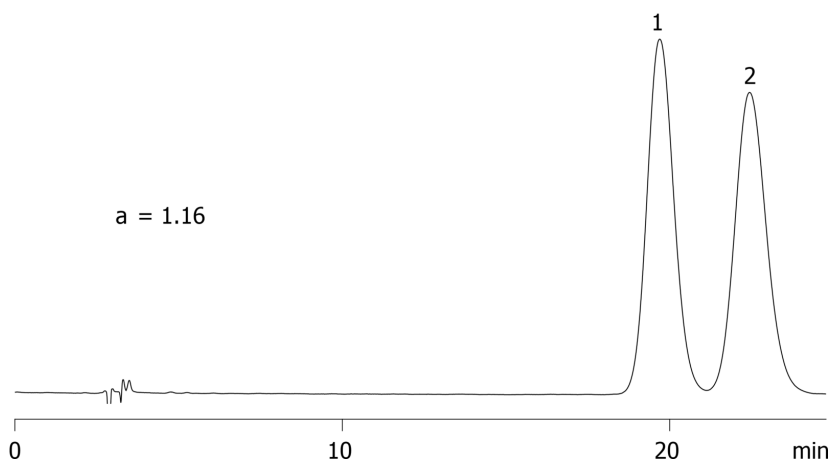

### ANALYTES:

- 1 Tropicamide  
Retention Time: 19.69 min
- 2 Tropicamide  
Retention Time: 22.44 min

