

Warfarin on Lux 5µm Cellulose-2 in NP

Column: Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4457-E0

Elution Type: Isocratic

Eluent A: Hexane / IPA (60:40) Formic acid (0.1%)

Gradient

Step No.	Time (min)	Pct A
1	0	100

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: Circular Dichroism @ 0.0000000000 (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, Lux Cellulose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8366

- Holder Part No.: KJ0-4282

20271



Products used in this application:

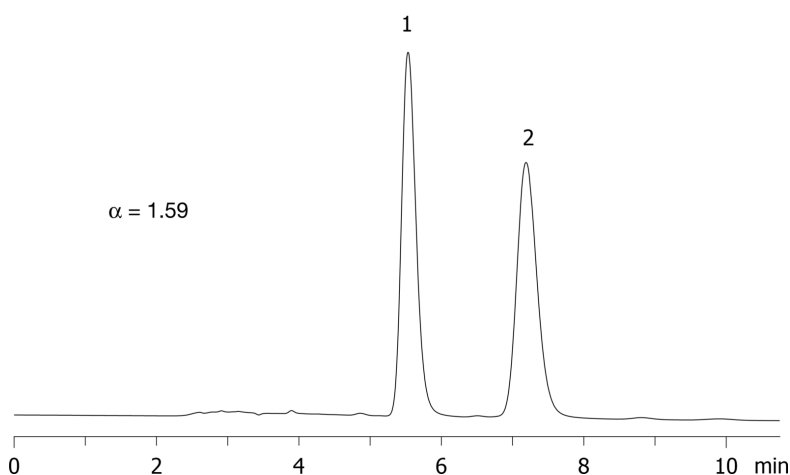

ANALYTES:

1 Warfarin

Retention Time: 5.53 min

2 Warfarin

Retention Time: 7.19 min

