

Ketamine on Lux 5µm Cellulose-3 in RP

Column: Lux® 5 µm Cellulose-3, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4493-E0

Elution Type: Isocratic

Eluent A: MeOH / 20mM NH₄HCO₃+0.1%DEA (90:10) None

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: Circular Dichroism @ 0.0000000000 (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-3 4 x 3.0 mm, 10/Pk Part No.: AJ0-8622


- Holder Part No.: KJ0-4282



Products used in this application:



20287

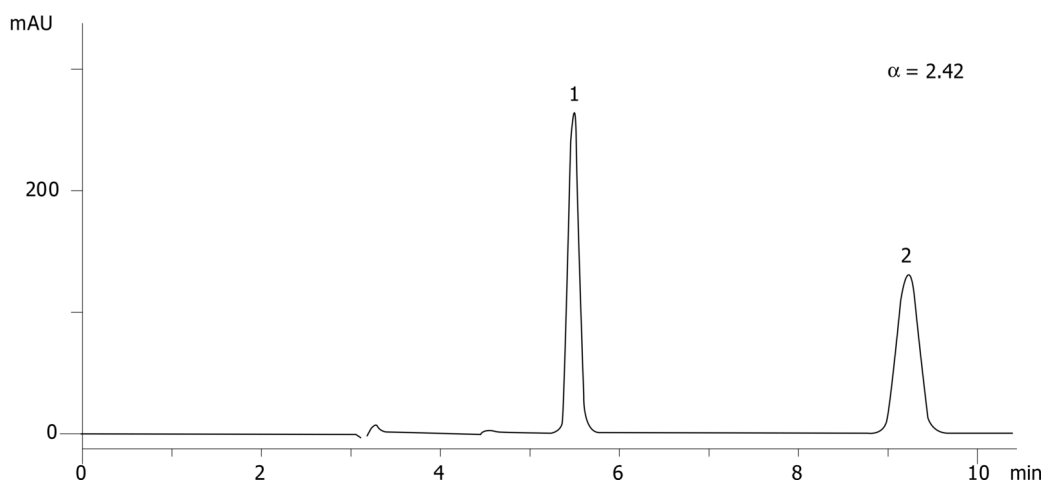

ANALYTES:

1 Ketamine

Retention Time: 5.49 min

2 Ketamine

Retention Time: 9.22 min

