

Warfarin on Lux 5µm Cellulose-3 in RP

Column: Lux® 5 µm Cellulose-3, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4493-E0

Elution Type: Isocratic

Eluent A: MeOH / 0.1% Formic Acid (90:10) None

Gradient

Step No. **Time (min)** **Pct A**

Profile: **1** **0** **100**

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: Circular Dichroism @ 0.0000000000 (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-3 4 x 3.0 mm, 10/Pk Part No.: AJ0-8622


- Holder Part No.: KJ0-4282



Products used in this application:



20295

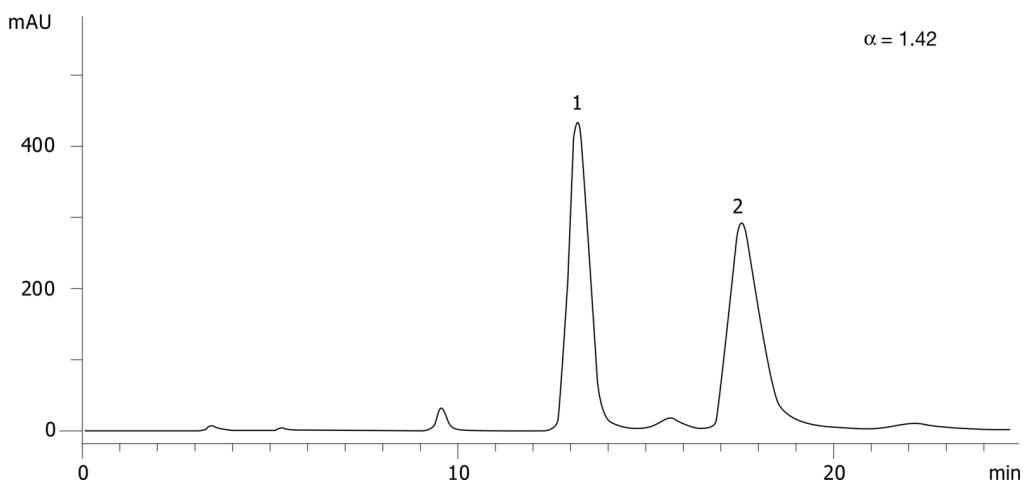

ANALYTES:

1 Warfarin

Retention Time: 13.18 min

2 Warfarin

Retention Time: 17.56 min

