

Diltiazem on Lux 5µm Amylose-2 in NP

Column: N/A

Dimensions: N/A

Order No: N/A

Elution Type: N/A

Eluent A: Hexane / EtOH (60:40) DEA (0.1%)

Gradient

Step No.	Time (min)	Pct A
1	0	100

Profile:

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: Circular Dichroism @ 0.0000000000 (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Amylose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8470

- Holder Part No.: KJ0-4282

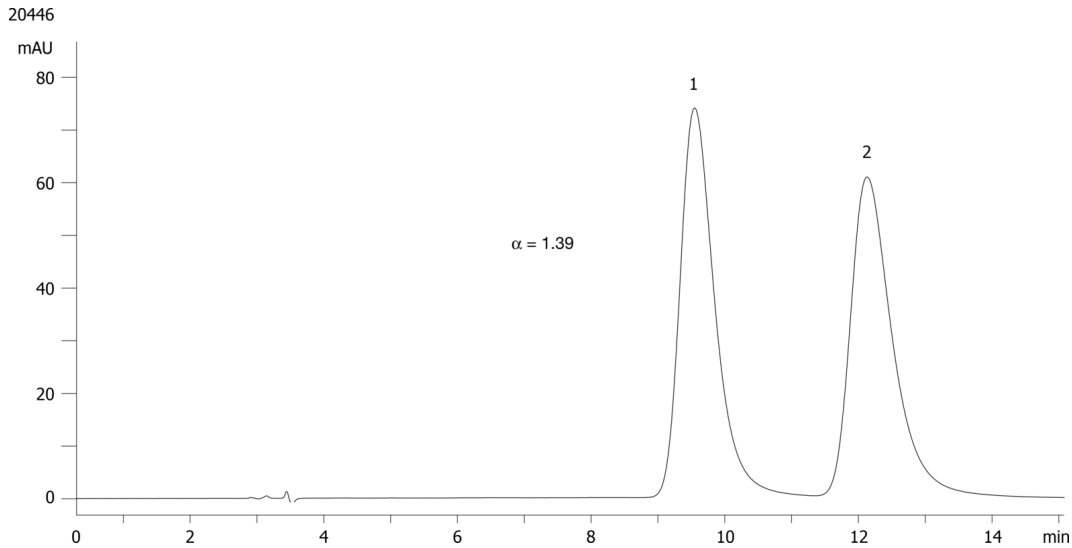

ANALYTES:

1 Diltiazem

Retention Time: 9.55 min

2 Diltiazem

Retention Time: 12.13 min

