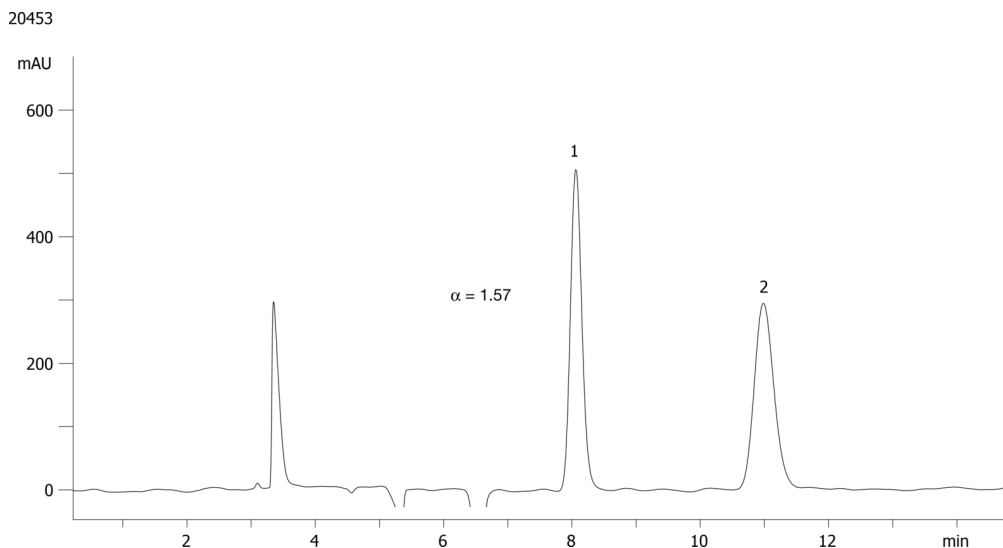


**Fenoprofen on Lux 5µm Amylose-2 in NP****Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Hexane / EtOH (95:5) Formic acid (0.1%)**Gradient**      **Step No.**      **Time (min)**      **Pct A****Profile:**              **1**              0              100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** Circular Dichroism @ 0.0000000000 (ambient)**Analyst Note:** SecurityGuard<sup>TM</sup> Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Amylose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8470

- Holder Part No.: KJ0-4282

**ANALYTES:****1** Fenoprofen

Retention Time: 8.06 min

**2** Fenoprofen

Retention Time: 10.98 min

