

Bambuterol on Lux 5µm Amylose-2 in NP

Column: N/A

Dimensions: N/A

Order No: N/A

Elution Type: N/A

Eluent A: Hexane / EtOH (60:40) DEA (0.1%)

Gradient **Step No.** **Time (min)** **Pct A**

Profile: **1** 0 100

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: Circular Dichroism @ 0.0000000000 (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Amylose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8470

- Holder Part No.: KJ0-4282

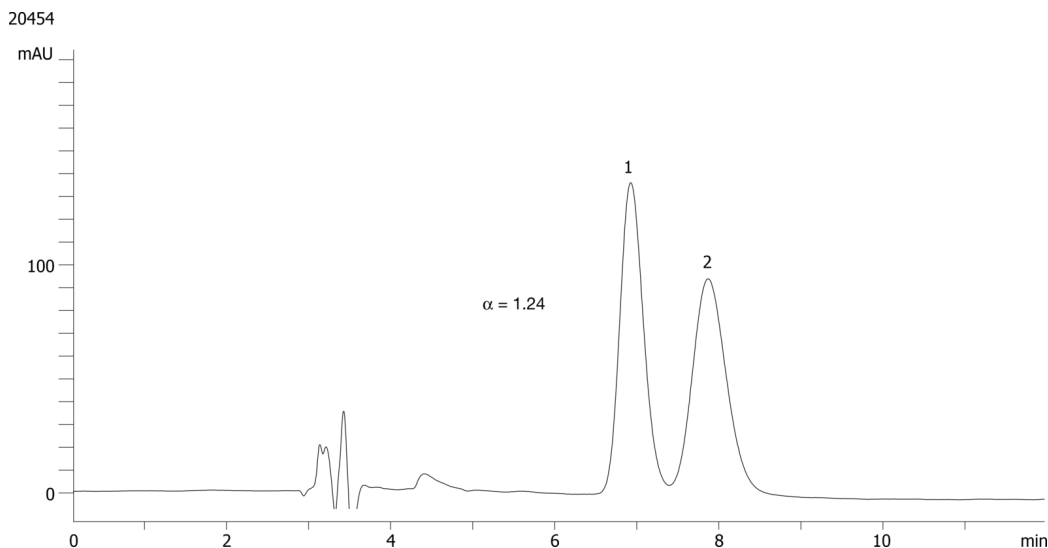

ANALYTES:

1 Bambuterol

Retention Time: 6.92 min

2 Bambuterol

Retention Time: 7.87 min

