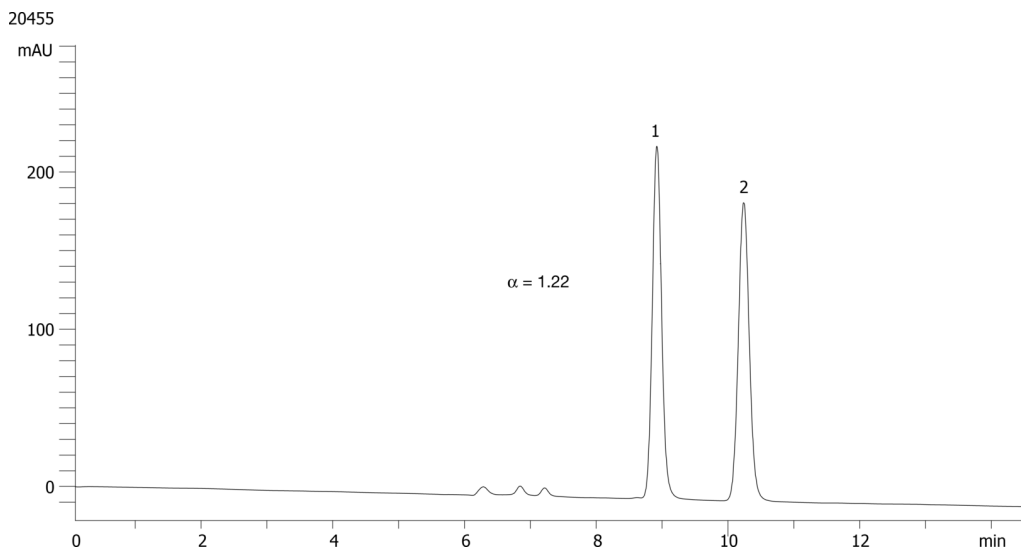


trans-Stilbene oxide on Lux 5 μ m Amylose-2 in PO**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** MeOH / IPA (90:10) DEA (0.1%)**Gradient** **Step No.** **Time (min)** **Pct A****Profile:** **1** 0 100**Flow Rate:** 0.5 mL/min**Col. Temp.:** ambient**Detection:** Circular Dichroism @ 0.0000000000 (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Amylose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8470

- Holder Part No.: KJ0-4282

**ANALYTES:****1** trans-Stilbene oxide

Retention Time: 8.92 min

2 trans-Stilbene oxide

Retention Time: 10.24 min

