

## Warfarin on Lux 5µm Amylose-2 in NP

**Column:** N/A

**Dimensions:** N/A

**Order No:** N/A

**Elution Type:** N/A

**Eluent A:** Hexane / EtOH (60:40) Formic acid (0.1%)

**Gradient**      **Step No.**      **Time (min)**      **Pct A**

**Profile:**              **1**              0              100

**Flow Rate:** 1 mL/min


**Col. Temp.:** ambient

**Detection:** Circular Dichroism @ 0.0000000000 (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Amylose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8470

- Holder Part No.: KJ0-4282

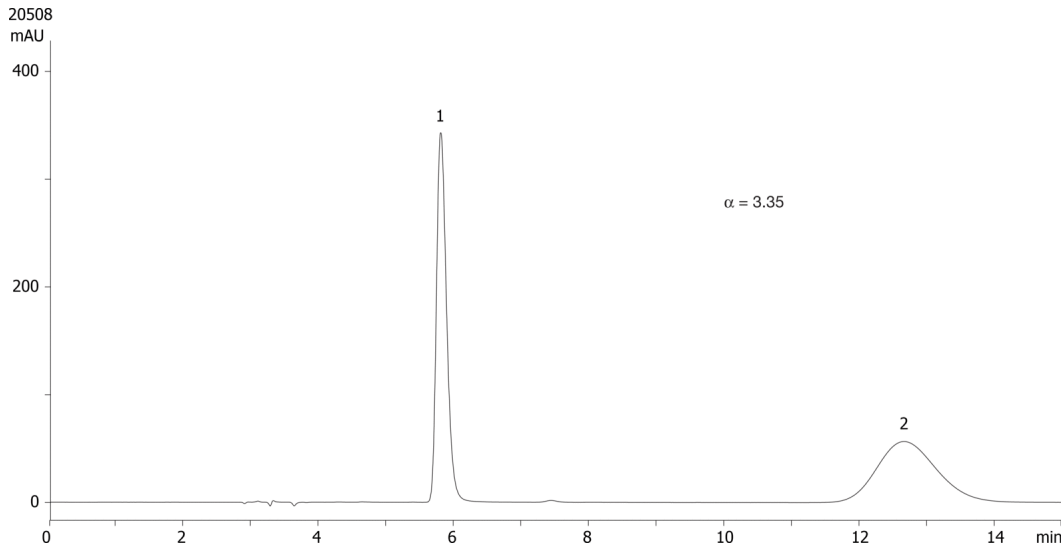

### ANALYTES:

**1** Warfarin

Retention Time: 5.81 min

**2** Warfarin

Retention Time: 12.66 min

