

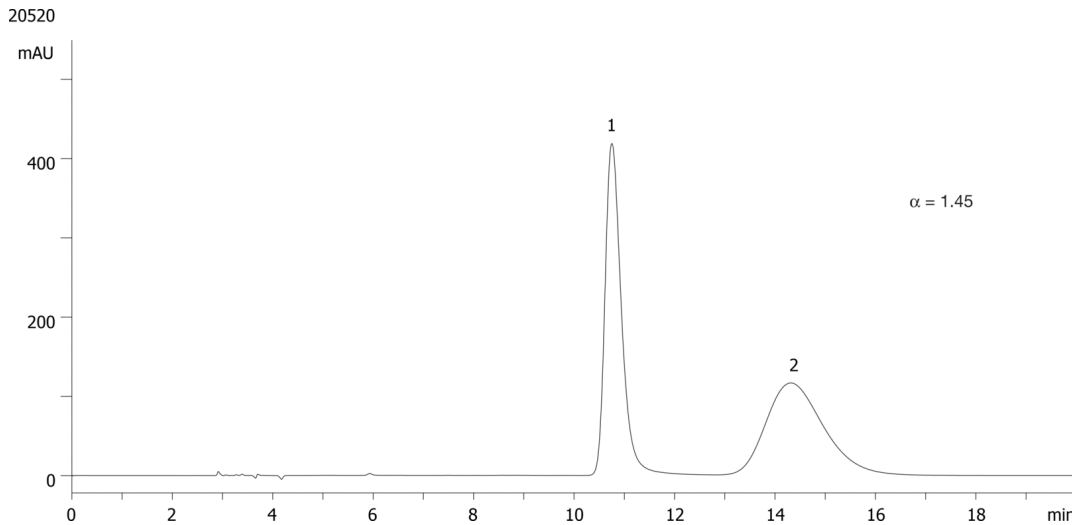
FMOC-Lys-(boc)-OH on Lux 5 μ m Amylose-2 in NP**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Hexane / IPA (80:20) Formic acid (0.1%)

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** Circular Dichroism @ 0.0000000000 (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Amylose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8470

- Holder Part No.: KJ0-4282

**ANALYTES:****1** FMOC-Lys-(Boc)-OH

Retention Time: 10.75 min

2 FMOC-Lys-(Boc)-OH

Retention Time: 14.31 min

