


Overlay: Yarra 3µm SEC-3000 versus Zenix SEC-300 Low Salt Buffer Competitive Comparisons

Column: Yarra™ 3 µm SEC-3000, LC Column 300 x 7.8 mm, Ea
Dimensions: 300 x 7.8 mm ID
Order No: 00H-4513-K0
Elution Type: Isocratic
Eluent A: 100 mM Sodium Phosphate in Water pH 6.8 with 0.025% NaN₃
Gradient

Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Circular Dichroism @ 0.0000000000 (ambient)

20612

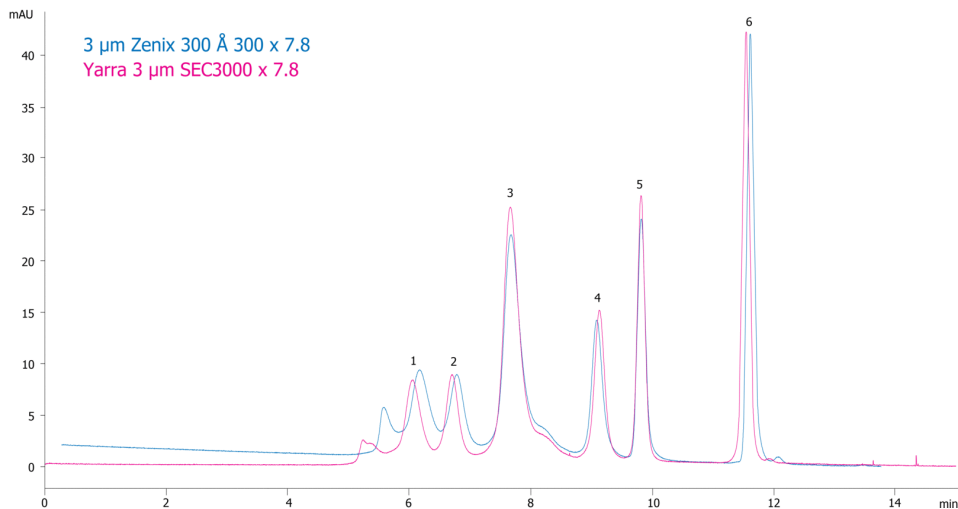

ANALYTES:

- 1** Thyroglobulin
Retention Time: 6 min
- 2** IgA
Retention Time: 6.7 min
- 3** IgG
Retention Time: 7.7 min
- 4** Ovalbumin
Retention Time: 9.2 min
- 5** Myoglobin
Retention Time: 9.8 min
- 6** Uridine
Retention Time: 11.5 min



Products used in this application:

