

Overlay: Yarra 3µm SEC-2000 versus Zenix SEC-150 low salt buffer PMQC Standard

Column: Yarra™ 3 µm SEC-2000, LC Column 300 x 4.6 mm, Ea

Dimensions: 300 x 4.6 mm ID

Order No: 00H-4512-E0

Elution Type: Isocratic


Eluent A: 100 mM Sodium Phosphate in Water pH 6.8 with 0.025% NaN₃

Gradient Profile:	Step No.	Time (min)	Pct A
	1	0	100

Flow Rate: 0.35 mL/min

Col. Temp.: ambient

Detection: Circular Dichroism @ 0.0000000000 (ambient)

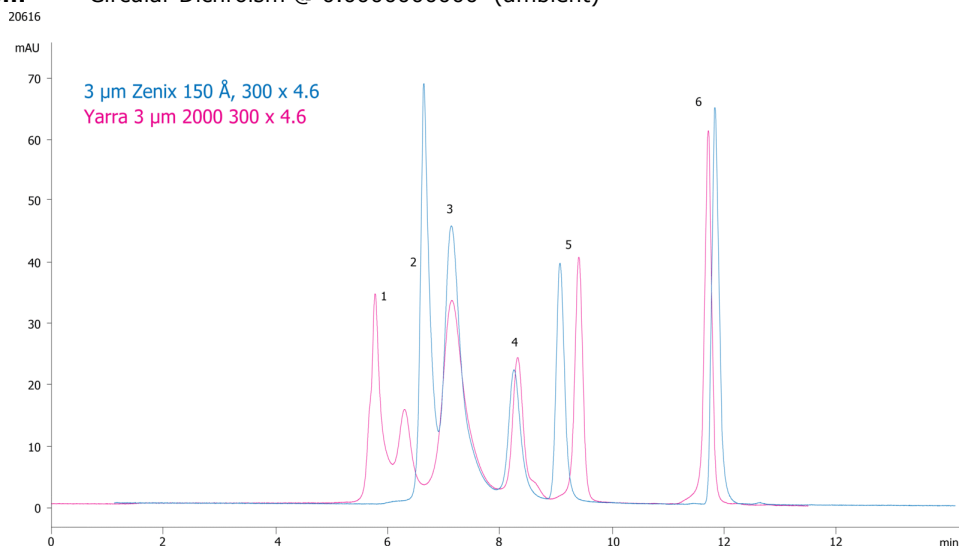

Products used in this application:



ANALYTES:

- 1** Thyroglobulin
Retention Time: 5.4 min
- 2** IgA
Retention Time: 5.9 min
- 3** IgG
Retention Time: 6.7 min
- 4** Ovalbumin
Retention Time: 8.2 min
- 5** Myoglobin
Retention Time: 9 min
- 6** Uridine
Retention Time: 11.3 min

