

Calibration curves for Yarra phases SEC-2000, SEC-3000, SEC-4000

Column: Yarra™ 3 µm SEC-2000, LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-4512-K0

Elution Type: Isocratic

Eluent A: 100 mM Sodium Chloride pH 6.8

Gradient Profile:

Step No.	Time (min)	Pct A
1	0	100


Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: Circular Dichroism @ 0.0000000000 (ambient)
20638



Products used in this application:

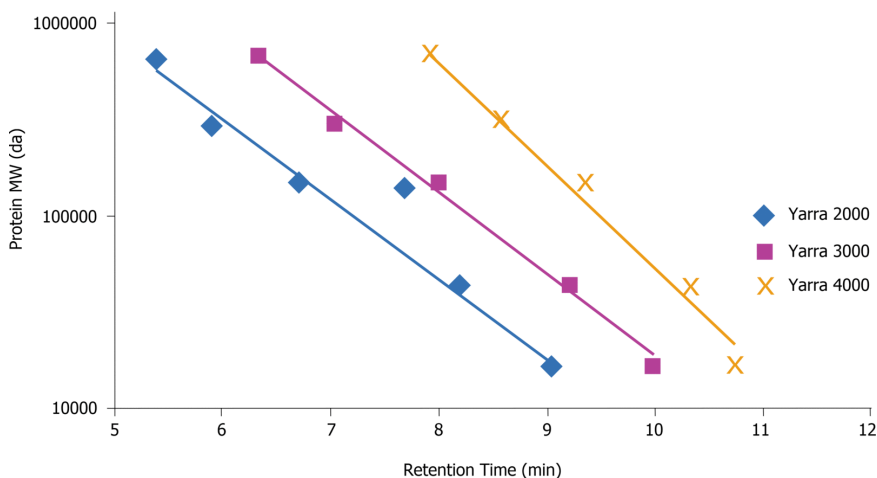

ANALYTES:

- 1 Thyroglobulin
- 2 IgA
- 3 IgG
- 4 Ovalbumin
- 5 Myoglobin
- 6 Uridine

