

Bath Salts analysis using Kinetex 2.6u C18 50x3.0mm

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 3 mm, Ea

Dimensions: 50 x 3 mm ID

Order No: 00B-4462-Y0

Elution Type: Gradient

Eluent A: 10mM Ammonium Formate

Eluent B: Methanol with 0.1% Formic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	3	0	100
	3	4	0	100
	4	4.1	95	5

Flow Rate: 500


Col. Temp.: ambient

Detection: Circular Dichroism @ 0.0000000000 (ambient)

20739



Products used in this application:

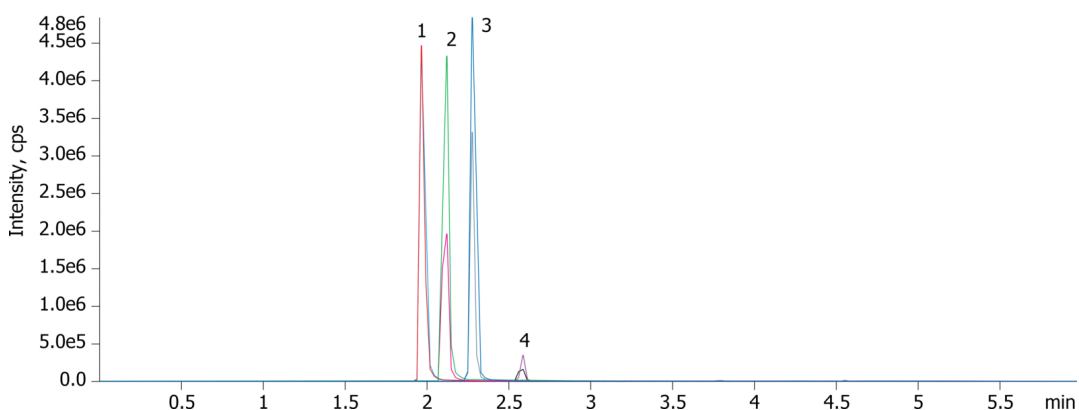

ANALYTES:

1 Methamphetamine

Retention Time: 1.12 min

2 Methylone

Retention Time: 1.97 min

3 Mephedrone

Retention Time: 2.28 min

4 MDPV

Retention Time: 2.69 min

