

Figure 3 - EuP Acetylsalicylic acid using Kinetex 2.6u C18 150x4.6mm

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4462-E0

Elution Type: Isocratic

Eluent A: 600/400/2 Water/Acetonitrile/H3PO4

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min

Col. Temp.: ambient


Detection: Circular Dichroism @ 0.0000000000 (ambient)



Products used in this application:



20768

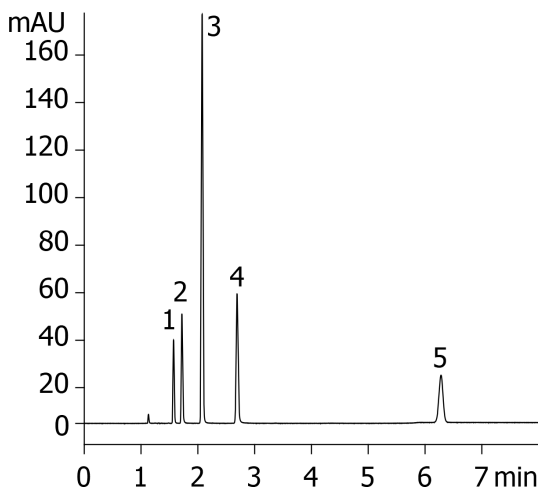

ANALYTES:

- 1 Impurity A
- 2 Impurity B
- 3 Impurity C
- 4 Acetylsalicylic acid
- 5 Impurity D

