

Pharmaceutical degradation using Kinetex 5u XB-C18 150x4.6mm

Column: Kinetex® 5µm XB-C18 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4605-E0

Elution Type: Gradient

Eluent A: Water with 0.1% TFA

Eluent B: Acetonitrile with 0.1% TFA

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	90	10
	2	20	30	70

Flow Rate: 1.2 mL/min


Col. Temp.: ambient

Detection: Circular Dichroism @ 0.0000000000 (ambient)

20784



Products used in this application:

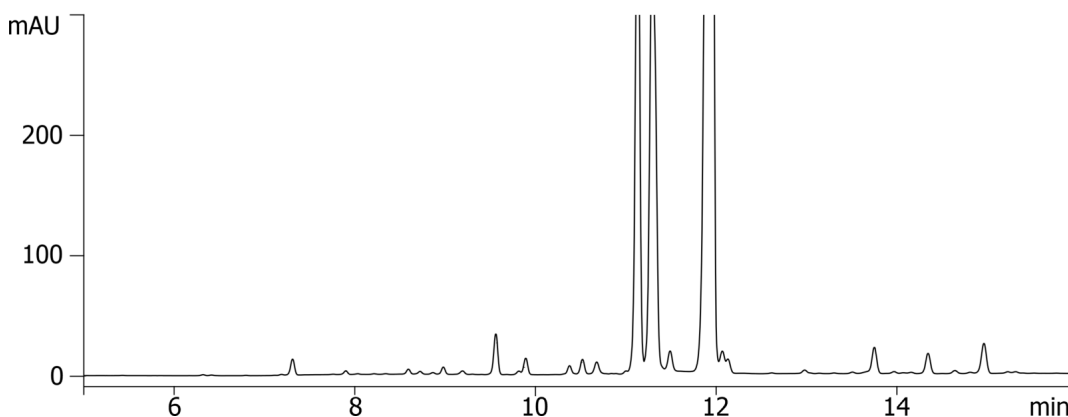

ANALYTES:

- 1 Mupirocin degradants

