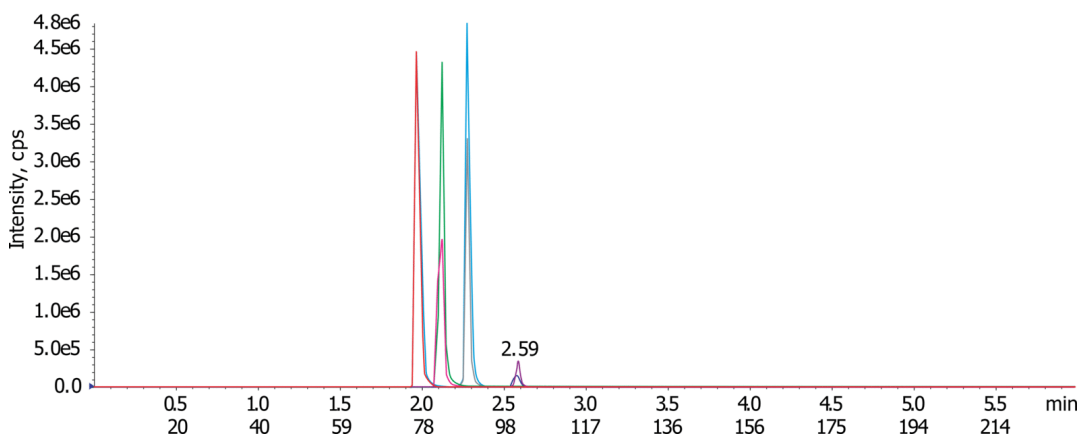


Bath salts analysis using Kinetex 2.6u C18 50x3.0mm**Column:** Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 3 mm, Ea**Dimensions:** 50 x 3 mm ID**Order No:** 00B-4462-Y0**Elution Type:** Gradient**Eluent A:** 10mM Ammonium formate**Eluent B:** Methanol with 0.1% Formic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	3	0	100
	3	4	0	100

Flow Rate: 500**Col. Temp.:** ambient**Detection:** Circular Dichroism @ 0.0000000000 (ambient)

20791

**Products used in this application:****ANALYTES:**

- 1 Methylone
Retention Time: 1.97 min
- 2 Methamphetamine
Retention Time: 2.12 min
- 3 Mephedrone
Retention Time: 2.28 min
- 4 MDPV
Retention Time: 2.69 min

