

Quechers(AOAC-LCMS) - Lettuce(iceberg) extract using Luna 3u C18(2) 150x3.0mm

Column: Luna® 3 µm C18(2) 100 Å, LC Column 150 x 3 mm, Ea

Dimensions: 150 x 3 mm ID

Order No: 00F-4251-Y0

Elution Type: Gradient

Eluent A: 5mM Formic Acid (no pH adjustment)

Eluent B: Methanol with 5mM Formic Acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	75	25
	2	15	10	90
	3	30	10	90
	4	30.01	75	25
	5	43	75	25

Flow Rate: 0.3 mL/min


Col. Temp.: ambient

Detection: Circular Dichroism @ 0.0000000000 (ambient)

20828



Products used in this application:

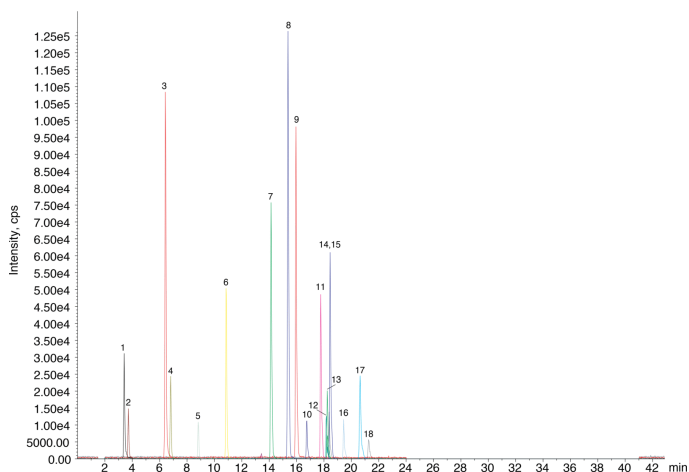

Quechers(AOAC-LCMS) - Lettuce(iceberg) extract using Luna 3u C18(2) 150x3.0mm

ANALYTES:

- 1** Pymetrozine
- 2** Methamidophos
- 3** Thiabendazole
- 4** Methomyl
- 5** Imidacloprid
- 6** Imazalil
- 7** Carbaryl
- 8** Atrazine
- 9** Azoxystrobin
- 10** Linuron
- 11** Cyprodinil
- 12** Parathion-d10
- 13** Kresoxim-methyl
- 14** Triphenyl phosphate
- 15** Tebuconazole
- 16** Chlorpyrifos-methyl
- 17** Ethion
- 18** Chlorpyrifos

