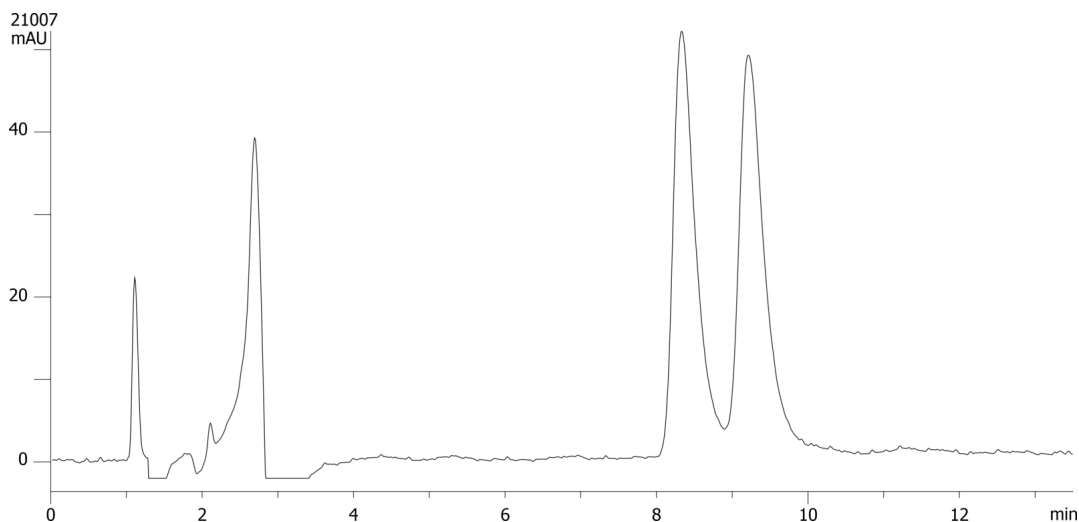


**Alprenolol on Lux 5µm Amylose-2 by SFC****Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** CO<sub>2</sub>/(Isopropanol + 0.5% IPA) 90:10 (v/v)**Gradient**      **Step No.**      **Time (min)**      **Pct A****Profile:**              **1**              0              100**Flow Rate:** 3 mL/min**Col. Temp.:** 30 °C**Detection:** Electrospray Mass Spec (ESMS) @ 0.0000000000 (30 °C)**Analyst Note:** Courtesy of Katrijn De Klerck, Debby Mangelings and Yvan Vander Heyden, Department of Analytical Chemistry and Pharmaceutical Technology, Vrije Universiteit Brussel, Brussels, Belgium.  
Results of this analysis yielded: Alpha >1.1 and Resolution >1.5SecurityGuard<sup>TM</sup> Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Amylose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8470

- Holder Part No.: KJ0-4282

**ANALYTES:****1** Alprenolol**2** Alprenolol