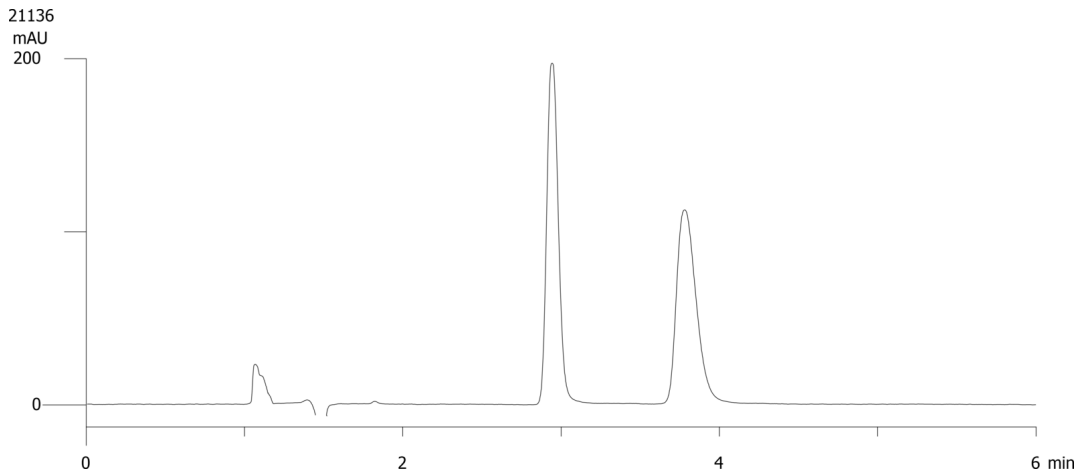


Flurbiprofen on Lux 5µm Amylose-2 by SFC**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** CO₂/(Methanol + 0.25% IPA + 0.25% TFA) 90:10 (v/v)**Gradient** **Step No.** **Time (min)** **Pct A****Profile:** **1** 0 100**Flow Rate:** 3 mL/min**Col. Temp.:** 30 °C**Detection:** Electrospray Mass Spec (ESMS) @ 0.0000000000 (30 °C)**Analyst Note:** Courtesy of Katrijn De Klerck, Debby Mangelings and Yvan Vander Heyden, Department of Analytical Chemistry and Pharmaceutical Technology, Vrije Universiteit Brussel, Brussels, Belgium.
Results of this analysis yielded: Alpha >1.1 and Resolution >1.5

SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Amylose-2, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8470

- Holder Part No.: KJ0-4282

**ANALYTES:****1** Flurbiprofen**2** Flurbiprofen