

Tetramisole on Lux 5µm Cellulose-4 by SFC**Column:** Lux® 5 µm Cellulose-4, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4491-E0**Elution Type:** Isocratic**Eluent A:** CO₂/(Methanol + 0.1% IPA + 0.1% TFA) 80:20 (v/v)

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 3 mL/min**Col. Temp.:** 30 °C**Detection:** Electrospray Mass Spec (ESMS) @ 0.0000000000 (30 °C)**Analyst Note:** Courtesy of Katrijn De Klerck, Debby Mangelings and Yvan Vander Heyden, Department of Analytical Chemistry and Pharmaceutical Technology, Vrije Universiteit Brussel, Brussels, Belgium.
Results of this analysis yielded: Alpha >1.1 and Resolution >1.5

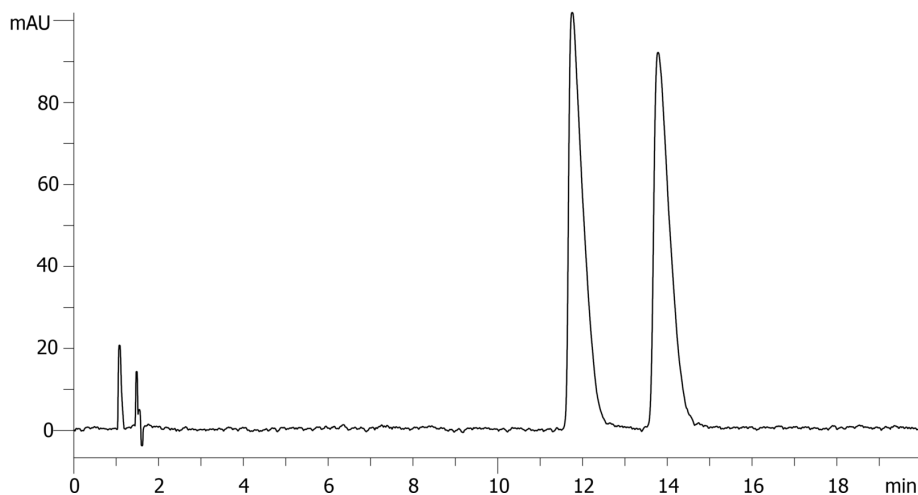
SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-4, 4 x 3.0mm ID, 10/Pk Part No.: AJ0-8627

- Holder Part No.: KJ0-4282

**Products used in this application:**

21406

**ANALYTES:****1** Tetramisole**2** Tetramisole