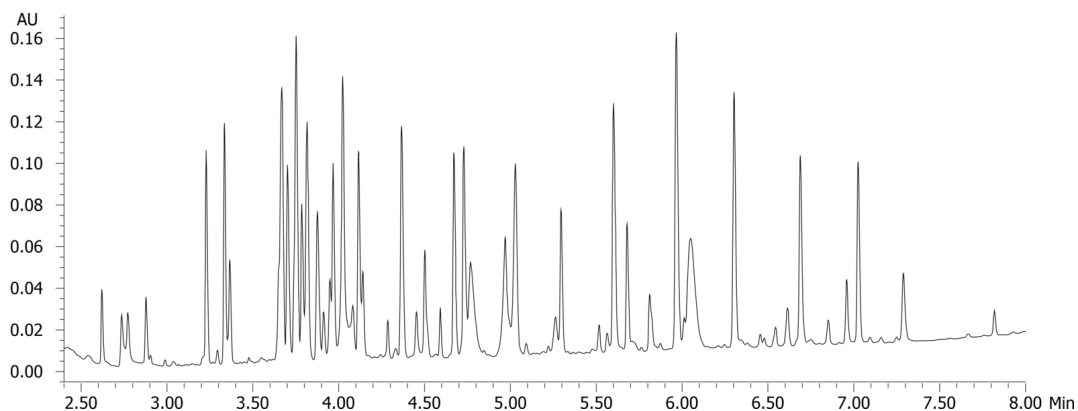


Rapid analysis of BSA digest using Kinetex 1.3u C18 50x2.1mm**Column:** Kinetex® 1.3 µm C18 100 Å, LC Column 50 x 2.1mm, Ea**Dimensions:** 50 x 2.1 mm ID**Order No:** 00B-4515-AN**Elution Type:** Gradient**Eluent A:** Water with 0.1% TFA**Eluent B:** Acetonitrile with 0.1% TFA

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	1	95	5
	3	11	5	95
	4	11.5	5	95
	5	11.6	95	5

Flow Rate: 0.5 mL/min**Col. Temp.:** 30 °C**Detection:** Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)

21459

**Products used in this application:****ANALYTES:**

- 1 Bovine Serum Albumin (BSA)

