

Overlay: Yarra 3µm SEC-2000 with different NaSO₄ MP


Column: Yarra™ 3 µm SEC-2000, LC Column 300 x 7.8 mm, Ea
Dimensions: 300 x 7.8 mm ID
Order No: 00H-4512-K0
Elution Type: Isocratic
Eluent A: 50 mM Sodium Phosphate in Water pH 6.8 + various salt concentrations
Gradient

Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)
21464



Products used in this application:

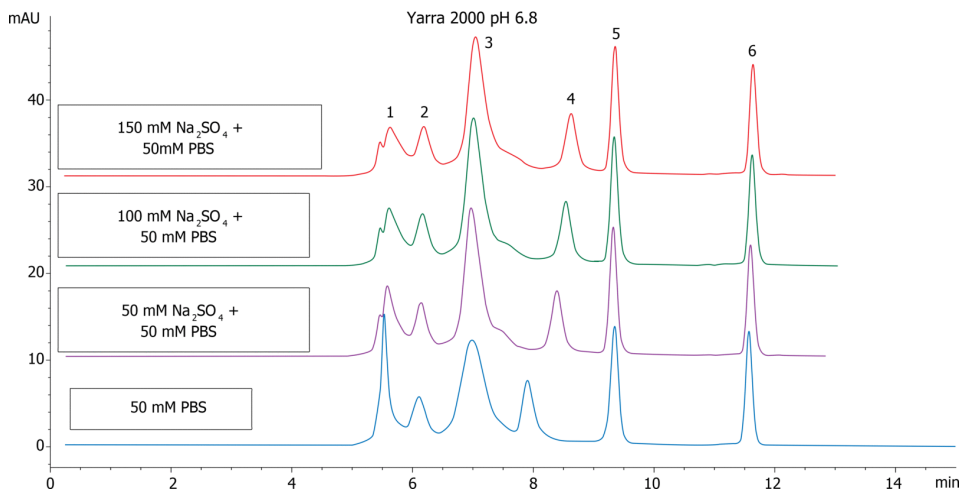

ANALYTES:

- 1** Thyroglobulin
Retention Time: 5.4 min
- 2** IgA
Retention Time: 5.9 min
- 3** IgG
Retention Time: 6.7 min
- 4** Ovalbumin
Retention Time: 8.2 min
- 5** Myoglobin
Retention Time: 9 min
- 6** Uridine
Retention Time: 11.3 min

